



SZABO SCANDIC

Part of Europa Biosite

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!
See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

mail@szabo-scandic.com

www.szabo-scandic.com

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic)



CD32B antibody [HMV315] HistoMAX™

Cat. No. GTX639939

Host	Rabbit
Clonality	Monoclonal
Isotype	IgG
Application	IHC-P
Reactivity	Human

Package

500 µl, 100 µl

PRODUCT

This antibody was validated on 76 different Normal Tissues by IHC-P.

Summary

[Go to Normal Tissue Gallery](#)[Go to Cancer Tissue Gallery](#)

APPLICATION

Application Note

*Optimal dilutions/concentrations should be determined by the researcher.

Suggested dilution	Recommended dilution
IHC-P	1:100-1:200

Note : Heat-induced antigen retrieval for 5 minutes in an autoclave at 121°C in pH 7.8 Target Retrieval Solution buffer.

Not tested in other applications.

Product Note

This antibody was raised against human CD32B Intracellular domain.

PROPERTIES

Form	Liquid
Buffer	PBS, 2% BSA
Preservative	0.025% ProClin 300
Storage	Store as concentrated solution. Centrifuge briefly prior to opening vial. For short-term storage (1-2 weeks), store at 4°C. For long-term storage, aliquot and store at -20°C or below. Avoid multiple freeze-thaw cycles.
Concentration	0.3 mg/ml
Immunogen	Synthetic peptide encompassing a sequence within the Intracellular domain of human CD32B. The exact sequence is proprietary.
Purification	Affinity purified by Protein A.
Conjugation	Unconjugated



For full product information, images and publications, please visit our [website](#).


Date 2024 / 12 / 12 Page 1 of 2

For laboratory research use only. Not for any clinical, therapeutic, or diagnostic use in humans or animals. Not for animal or human consumption.

Note

Purchasers shall not, and agree not to enable third parties to, analyze, copy, reverse engineer or otherwise attempt to determine the structure or sequence of the product.


DATA IMAGES

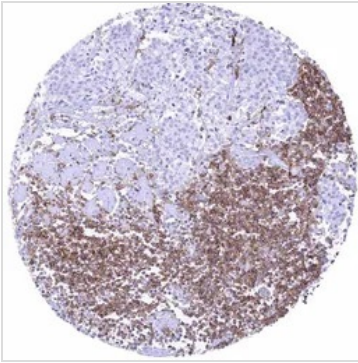

GTX639939 IHC-P Image

IHC-P analysis of human prostatic adenocarcinoma (PAC) tissue with (Gleason 5+5=10) tissue using GTX639939 CD32B antibody [HMV315] HistoMAX™.

CD32b negative adenocarcinoma (Gleason 5+5=10). CD32b staining occurs in tumor associated B lymphocytes.

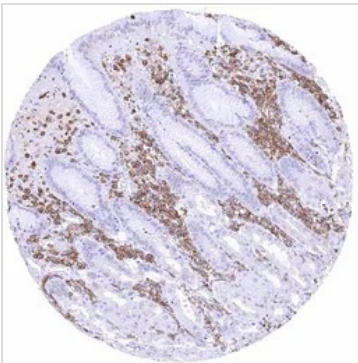

GTX639939 IHC-P Image

IHC-P analysis of human stomach tissue using GTX639939 CD32B antibody [HMV315] HistoMAX™.

Strong CD32B staining of many inflammatory cells.



GTX639939 IHC-P Image

IHC-P analysis of human diffuse large B-cell lymphoma (DLBCL) from lymph node tissue using GTX639939 CD32B antibody [HMV315] HistoMAX™.

Diffuse large B-cell lymphoma with strong CD32B staining of tumor cells.



For full product information, images and publications, please visit our [website](https://www.genetex.com).