



SZABO SCANDIC

Part of Europa Biosite

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!
See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

mail@szabo-scandic.com

www.szabo-scandic.com

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 

TIGIT antibody [HMV322] HistoMAX™

Cat. No. GTX639946

Host	Rabbit
Clonality	Monoclonal
Isotype	IgG
Application	IHC-P
Reactivity	Human

Package

500 µl, 100 µl

PRODUCT

This antibody was validated on 76 different Normal Tissues by IHC-P.

Summary

[Go to Normal Tissue Gallery](#)[Go to Cancer Tissue Gallery](#)

APPLICATION

Application Note

*Optimal dilutions/concentrations should be determined by the researcher.

Suggested dilution	Recommended dilution
IHC-P	1:100-1:200

Note : Heat-induced antigen retrieval for 5 minutes in an autoclave at 121°C in pH 7.8 Target Retrieval Solution buffer.

Not tested in other applications.

Product Note

This antibody was raised against human TIGIT Intracellular domain.

PROPERTIES

Form	Liquid
Buffer	PBS, 2% BSA
Preservative	0.025% ProClin 300
Storage	Store as concentrated solution. Centrifuge briefly prior to opening vial. For short-term storage (1-2 weeks), store at 4°C. For long-term storage, aliquot and store at -20°C or below. Avoid multiple freeze-thaw cycles.
Concentration	0.13 mg/ml
Immunogen	Recombinant protein encompassing a sequence within the Intracellular domain of human TIGIT. The exact sequence is proprietary.
Purification	Affinity purified by Protein A.
Conjugation	Unconjugated

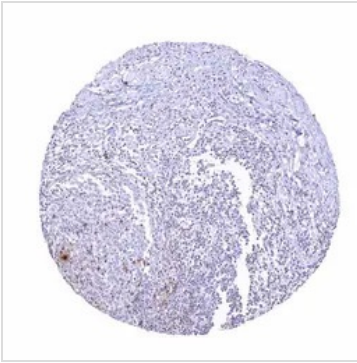


For full product information, images and publications, please visit our [website](#).

For laboratory research use only. Not for any clinical, therapeutic, or diagnostic use in humans or animals. Not for animal or human consumption.

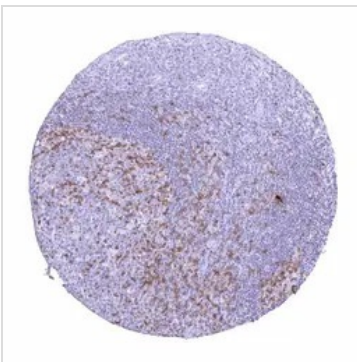
Note

Purchasers shall not, and agree not to enable third parties to, analyze, copy, reverse engineer or otherwise attempt to determine the structure or sequence of the product.

DATA IMAGES

GTX639946 IHC-P Image

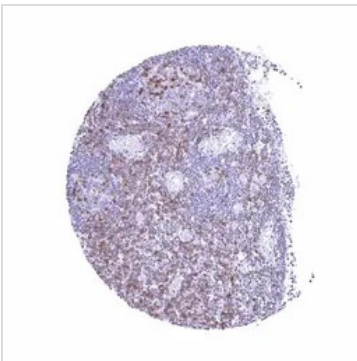
IHC-P analysis of human mucosa from appendix tissue using GTX639946 TIGIT antibody [HMOV322] HistoMAX™.

Distinct TIGIT staining of a subset of lymphocytes.


GTX639946 IHC-P Image

IHC-P analysis of human tonsil tissue using GTX639946 TIGIT antibody [HMOV322] HistoMAX™.

Most TIGIT positive lymphocytes are in the interfollicular area while the strongest TIGIT staining occurs in a subset of lymphocytes in the germinal centre.


GTX639946 IHC-P Image

IHC-P analysis of human lymph node tissue using GTX639946 TIGIT antibody [HMOV322] HistoMAX™.

Most TIGIT positive lymphocytes are in the interfollicular area while the strongest TIGIT staining occurs in the few labeled lymphocytes in a small germinal centre.



For full product information, images and publications, please visit our [website](#).