



SZABO SCANDIC

Part of Europa Biosite

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!
See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

mail@szabo-scandic.com

www.szabo-scandic.com

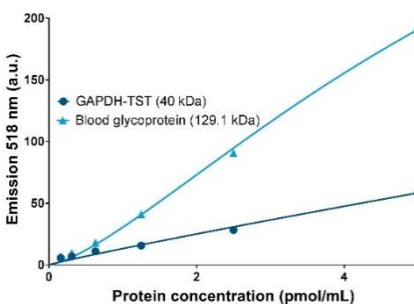
[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 

Data Sheet

Strep-Tactin®XT DY-488

Cat. No.: 2-1562-050

Version: 5.0
Revision Date: 11.01.2024

Description	Strep-Tactin®XT DY-488 is a streptavidin variant conjugated with the fluorescent dye DY-488. It can be used for detection of N-, C-terminally as well as internally tagged Strep-tag®II or Twin-Strep-tag® proteins in applications like ELISA, FACS, immunohistochemistry, and immunocytochemistry.																						
Specificity	Strep-tag®II, Twin-Strep-tag®, cross reactivity: biotin																						
Product Type	Recombinant protein																						
Form	Suspension in PBS; contains 0.01% NaN ₃ , 1 mg/ml BSA																						
Concentration	0.5 mg/ml																						
Amount	100 µl																						
Stability	6 months after shipping (If the products contain a precipitate, we recommend centrifugation before use.)																						
Storage	2-8 °C, protect from light																						
Shipping	Room temperature																						
Hazards	Product is not classified as hazardous according to (EC) No 1272/2008 [CLP]. A Material Safety Data Sheet is provided.																						
Absorption/Emission	493/518 nm																						
Application	<div>  <table border="1"> <caption>Approximate data points from the graph</caption> <thead> <tr> <th>Protein concentration (pmol/mL)</th> <th>GAPDH-TST (40 kDa) Emission 518 nm (a.u.)</th> <th>Blood glycoprotein (129.1 kDa) Emission 518 nm (a.u.)</th> </tr> </thead> <tbody> <tr> <td>0</td> <td>0</td> <td>0</td> </tr> <tr> <td>1</td> <td>~10</td> <td>~40</td> </tr> <tr> <td>2</td> <td>~20</td> <td>~80</td> </tr> <tr> <td>3</td> <td>~30</td> <td>~120</td> </tr> <tr> <td>4</td> <td>~40</td> <td>~160</td> </tr> <tr> <td>5</td> <td>~50</td> <td>~200</td> </tr> </tbody> </table> </div> <div> <p>FACS/Flow cytometry: Direct: dilute 1:5000 Indirect (via tagged protein): 75 ng (1.5 µl of 1:10 dilution) conjugate per 200 ng of a 50 kDa protein/50 ng of a 16 kDa protein for 5 x 10⁶ total cells</p> <p>ELISA: use a black plate with clear bottom, dilute between 1:500 and 1:2 000; start with 1:1 000.</p> <p>The detection limit may vary depending on protein size. The data provided on the left was obtained with Twin-Strep-tag® fused proteins in an ELISA, with a Strep-Tactin®XT DY-488 dilution of 1:1 000.</p> </div>		Protein concentration (pmol/mL)	GAPDH-TST (40 kDa) Emission 518 nm (a.u.)	Blood glycoprotein (129.1 kDa) Emission 518 nm (a.u.)	0	0	0	1	~10	~40	2	~20	~80	3	~30	~120	4	~40	~160	5	~50	~200
Protein concentration (pmol/mL)	GAPDH-TST (40 kDa) Emission 518 nm (a.u.)	Blood glycoprotein (129.1 kDa) Emission 518 nm (a.u.)																					
0	0	0																					
1	~10	~40																					
2	~20	~80																					
3	~30	~120																					
4	~40	~160																					
5	~50	~200																					

For research use only

Trademark information

The owners of trademarks marked by "®" or "TM" are identified at <http://www.iba-lifesciences.com/patents.html>. Registered names, trademarks, etc. used in this document, even when not specifically marked as such, are not to be considered unprotected by law.

Important licensing information

This product is covered by intellectual property (IP) rights and on completion of the sale IBA Lifesciences grants respective Limited Use Label Licenses to purchaser. IP rights and Limited Use Label Licenses for said technology are further described and identified at <http://www.iba-lifesciences.com/patents.html> or upon inquiry at info@iba-lifesciences.com or at IBA Lifesciences GmbH, Rudolf-Wissell-Str. 28, 37079 Goettingen, Germany. By use of this product the purchaser accepts the terms and conditions of all applicable Limited Use Label Licenses.