



SZABO SCANDIC

Part of Europa Biosite

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!
See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC Handels GmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

mail@szabo-scandic.com

www.szabo-scandic.com

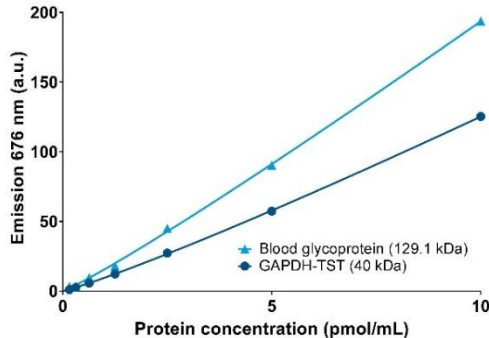
[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 

Data Sheet

Strep-Tactin®XT DY-649

Cat. No.: 2-1568-050

Version: 5.0
Revision Date: 11.01.2024

Description	Strep-Tactin®XT DY-649 is a streptavidin variant conjugated with the fluorescent dye DY-649. It can be used for detection of N-, C-terminally as well as internally tagged Strep-tag®II or Twin-Strep-tag® proteins in applications like ELISA, FACS, immunohistochemistry, and immunocytochemistry.																
Specificity	Strep-tag®II, Twin-Strep-tag®, cross reactivity: biotin																
Product Type	Recombinant protein																
Form	Suspension in PBS; contains 0.01% NaN ₃ , 1 mg/ml BSA																
Concentration	0.5 mg/ml																
Amount	100 µl																
Stability	6 months after shipping (If the products contain a precipitate, we recommend centrifugation before use.)																
Storage	2-8 °C, protect from light																
Shipping	Room temperature																
Hazards	Product is not classified as hazardous according to (EC) No 1272/2008 [CLP]. A Material Safety Data Sheet is provided.																
Absorption/Emission	654/672 nm																
Application	<div>  <table border="1"> <caption>Approximate data points from the graph</caption> <thead> <tr> <th>Protein concentration (pmol/mL)</th> <th>Blood glycoprotein (129.1 kDa) Emission (a.u.)</th> <th>GAPDH-TST (40 kDa) Emission (a.u.)</th> </tr> </thead> <tbody> <tr> <td>0</td> <td>0</td> <td>0</td> </tr> <tr> <td>2.5</td> <td>~45</td> <td>~25</td> </tr> <tr> <td>5</td> <td>~90</td> <td>~50</td> </tr> <tr> <td>10</td> <td>~180</td> <td>~100</td> </tr> </tbody> </table> </div> <div> <p>FACS/Flow cytometry: Direct: dilute 1:5000 Indirect (via tagged protein): 75 ng (1.5 µl auf a 1:10 dilution) conjugate per 200 ng of a 50 kDa protein/50 ng of a 16 kDa protein for 5 x 10⁶ total cells</p> <p>ELISA: use a black plate with clear bottom, dilute between 1:500 and 1:2 000; start with 1:500. The detection limit may vary depending on protein size. The data provided on the left was obtained with Twin-Strep-tag® fused proteins in an ELISA, with a Strep-Tactin®XT DY-649 dilution of 1:500.</p> </div>		Protein concentration (pmol/mL)	Blood glycoprotein (129.1 kDa) Emission (a.u.)	GAPDH-TST (40 kDa) Emission (a.u.)	0	0	0	2.5	~45	~25	5	~90	~50	10	~180	~100
Protein concentration (pmol/mL)	Blood glycoprotein (129.1 kDa) Emission (a.u.)	GAPDH-TST (40 kDa) Emission (a.u.)															
0	0	0															
2.5	~45	~25															
5	~90	~50															
10	~180	~100															

For research use only

Trademark information

The owners of trademarks marked by "®" or "TM" are identified at <http://www.iba-lifesciences.com/patents.html>. Registered names, trademarks, etc. used in this document, even when not specifically marked as such, are not to be considered unprotected by law.

Important licensing information

This product is covered by intellectual property (IP) rights and on completion of the sale IBA Lifesciences grants respective Limited Use Label Licenses to purchaser. IP rights and Limited Use Label Licenses for said technology are further described and identified at <http://www.iba-lifesciences.com/patents.html> or upon inquiry at info@iba-lifesciences.com or at IBA Lifesciences GmbH, Rudolf-Wissell-Str. 28, 37079 Goettingen, Germany. By use of this product the purchaser accepts the terms and conditions of all applicable Limited Use Label Licenses.