



# SZABO SCANDIC

Part of Europa Biosite

## Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!  
See the following pages for more information!



### Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

### Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

### SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

[mail@szabo-scandic.com](mailto:mail@szabo-scandic.com)

[www.szabo-scandic.com](http://www.szabo-scandic.com)

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 



IBA Lifesciences GmbH  
Rudolf-Wissell-Str. 28  
37079 Goettingen  
Germany  
Tel.: +49 (0) 551-5 06 72-0  
E-mail: [info@iba-lifesciences.com](mailto:info@iba-lifesciences.com)  
[www.iba-lifesciences.com](http://www.iba-lifesciences.com)

## Data Sheet

### pASG-IBA102

Cat. No.: 5-4102-001

Version: 3.0  
Revision Date: 26.07.2021

<b>Description</b>	<b>StarGate Acceptor Vector for bacterial expression. The expression cassette is under transcriptional control of the tetracycline promoter/operator. The expressed recombinant protein will be secreted into the periplasm.</b>
<b>Affinity tag</b>	Twin-Strep-tag® is fused to the C-terminus of the recombinant protein.
<b>Secretion</b>	The <i>ompA</i> signal sequence directs the expressed protein into the periplasmic space and will be cleaved off during the translocation process
<b>Cloning Strategy</b>	Cloning into StarGate Acceptor Vectors has to be done with the restriction enzyme Esp3I. There is no Multiple Cloning Site (MCS) available that can be used for the integration of the gene of interest instead (see manual).
<b>Expression strain</b>	Any <i>E. coli</i> strain. The <i>tet</i> -promoter works independently from the genetic background of <i>E. coli</i> .
<b>Bacterial Expression</b>	Expression is induced upon addition of 200 µg anhydrotetracycline per 1 liter <i>E. coli</i> shaking culture ( $A_{550} = 0.5$ ).
<b>Resistance</b>	Ampicillin
<b>Form</b>	5 µg, dissolved in 20 µl TE buffer, pH 8.0: 10 mM Tris/HCl, 1 mM EDTA
<b>Concentration</b>	250 ng/µl
<b>Stability</b>	12 months after shipping
<b>Storage</b>	recommended: 2-8 °C for frequent usage, -20 °C for long-term storage
<b>Shipping</b>	room temperature
<b>Hazards</b>	Product is not classified as hazardous according to (EC) No 1272/2008 [CLP]. A Material Safety Data Sheet is provided.

#### For research use only


##### Trademark information

The owners of trademarks marked by "®" or "TM" are identified at <http://www.iba-lifesciences.com/patents.html>. Registered names, trademarks, etc. used in this document, even when not specifically marked as such, are not to be considered unprotected by law.

##### Important licensing information


This product is covered by intellectual property (IP) rights and on completion of the sale IBA Lifesciences grants respective Limited Use Label Licenses to purchaser. IP rights and Limited Use Label Licenses for said technology are further described and identified at <http://www.iba-lifesciences.com/patents.html> or upon inquiry at [info@iba-lifesciences.com](mailto:info@iba-lifesciences.com) or at IBA Lifesciences GmbH, Rudolf-Wissell-Str. 28, 37079 Goettingen, Germany. By use of this product the purchaser accepts the terms and conditions of all applicable Limited Use Label Licenses.

Expression cassette of pASG-IBA102

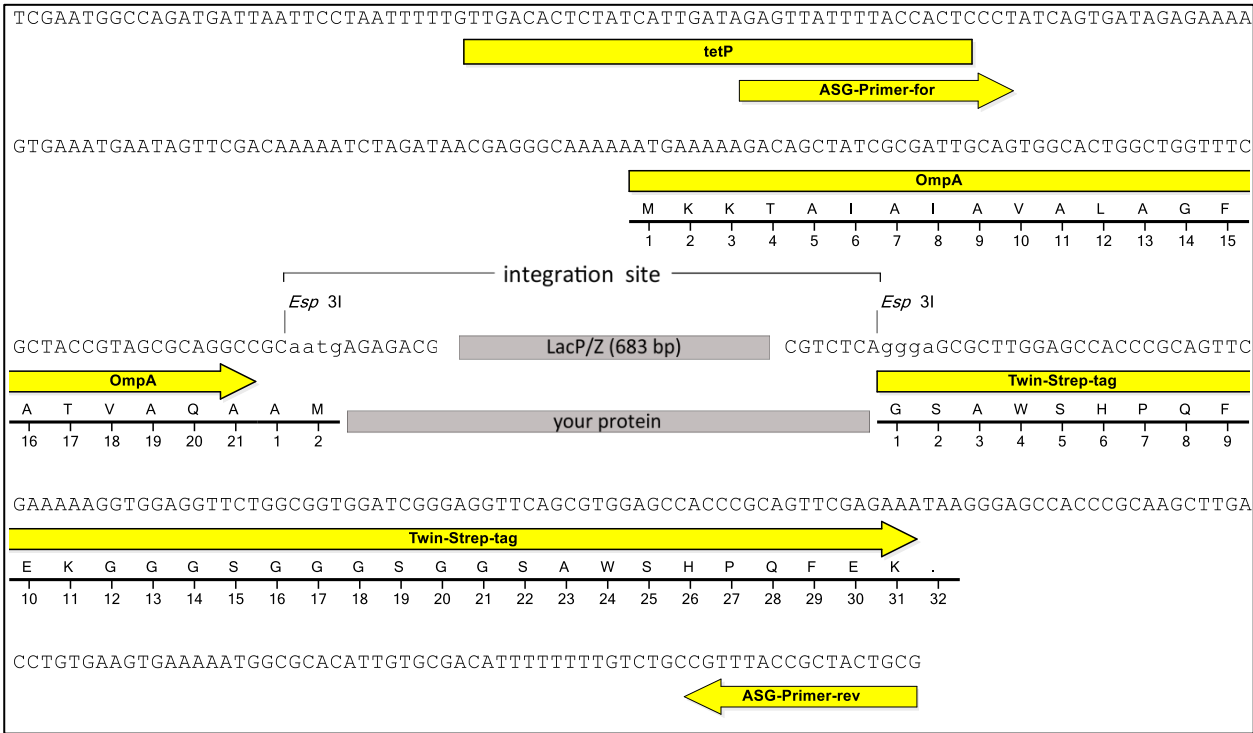

LacP/Z cassette = contains LacZ alpha fragment under control of a separate promoter, which allows alpha complementation of *LacZ* mutations such as *LacZΔM15* as in *E. coli* DH5α or TOP10.

your protein = after StarGate cloning using *Esp3I* your gene of interest will be located here

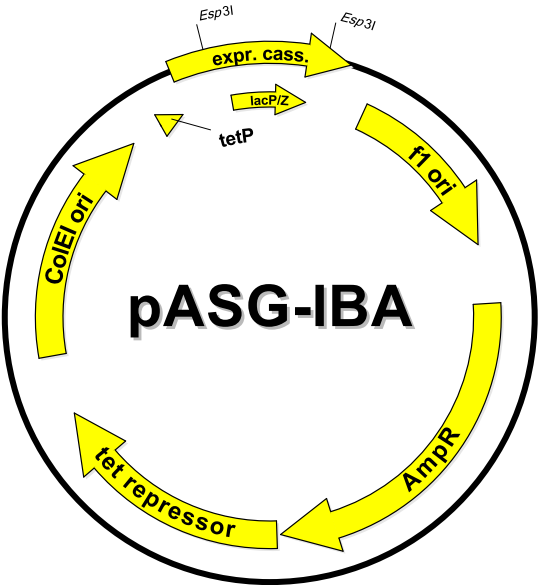

Features	from bp	to bp	Sequencing primer
f1 origin	13	451	ASG-Primer-for 5' - GAGTTATTTTACCACTCCCT -3'
AmpR resistance gene	600	1460	
Tet-repressor	1470	2093	
ColEI ori	2246	2834	
Tet promoter	2939	2975	ASG-Primer-rev 5' - CGCAGTAGCGGTAAACG -3'
forward primer binding site	2959	2978	
OmpA signal sequence	3041	3103	
LacZ alpha fragment	3335	3736	
Twin-Strep-tag®	3800	3892	
reverse primer binding site	3966	3982	
total vector length		3982	