



SZABO SCANDIC

Part of Europa Biosite

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!
See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

mail@szabo-scandic.com

www.szabo-scandic.com

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 



IBA Lifesciences GmbH
Rudolf-Wissell-Str. 28
37079 Goettingen
Germany
Tel.: +49 (0) 551-5 06 72-0
E-mail: info@iba-lifesciences.com
www.iba-lifesciences.com

Data Sheet

pESG-IBA105

Cat. No.: 5-4505-001

Version: 3.0
Revision Date: 27.07.2021

Description	StarGate Acceptor Vector for stable and non-replicative transient expression in mammalian cells. It carries the human cytomegalovirus (CMV) immediate-early promoter for high-level expression and the SV40 ori for episomal replication in cell lines that are latently infected with SV40 or that express the SV40 large T antigen (e.g., HEK293T, COS-1, COS-7). The expressed recombinant protein will be localized in the cytoplasm.
Affinity tag	Twin-Strep-tag® is fused to the N-terminus of the recombinant protein.
Cloning Strategy	Cloning into StarGate Acceptor Vectors has to be done with the restriction enzyme Esp3I. There is no Multiple Cloning Site (MCS) available that can be used for the integration of the gene of interest instead (see manual).
Resistance	Ampicillin: for selection of transformed <i>E. coli</i> cells Neomycin: for selection of stable cell lines
Form	5 µg, dissolved in 20 µl TE buffer, pH 8.0: 10 mM Tris/HCl, 1 mM EDTA
Concentration	250 ng/µl
Stability	12 months after shipping
Storage	recommended: 2-8 °C for frequent usage, -20 °C for long-term storage
Shipping	room temperature
Hazards	Product is not classified as hazardous according to (EC) No 1272/2008 [CLP]. A Material Safety Data Sheet is provided.

For research use only


Trademark information

The owners of trademarks marked by “®” or “TM” are identified at <http://www.iba-lifesciences.com/patents.html>. Registered names, trademarks, etc. used in this document, even when not specifically marked as such, are not to be considered unprotected by law.

Important licensing information


This product is covered by intellectual property (IP) rights and on completion of the sale IBA Lifesciences grants respective Limited Use Label Licenses to purchaser. IP rights and Limited Use Label Licenses for said technology are further described and identified at <http://www.iba-lifesciences.com/patents.html> or upon inquiry at info@iba-lifesciences.com or at IBA Lifesciences GmbH, Rudolf-Wissell-Str. 28, 37079 Goettingen, Germany. By use of this product the purchaser accepts the terms and conditions of all applicable Limited Use Label Licenses.

Expression cassette of pESG-IBA105

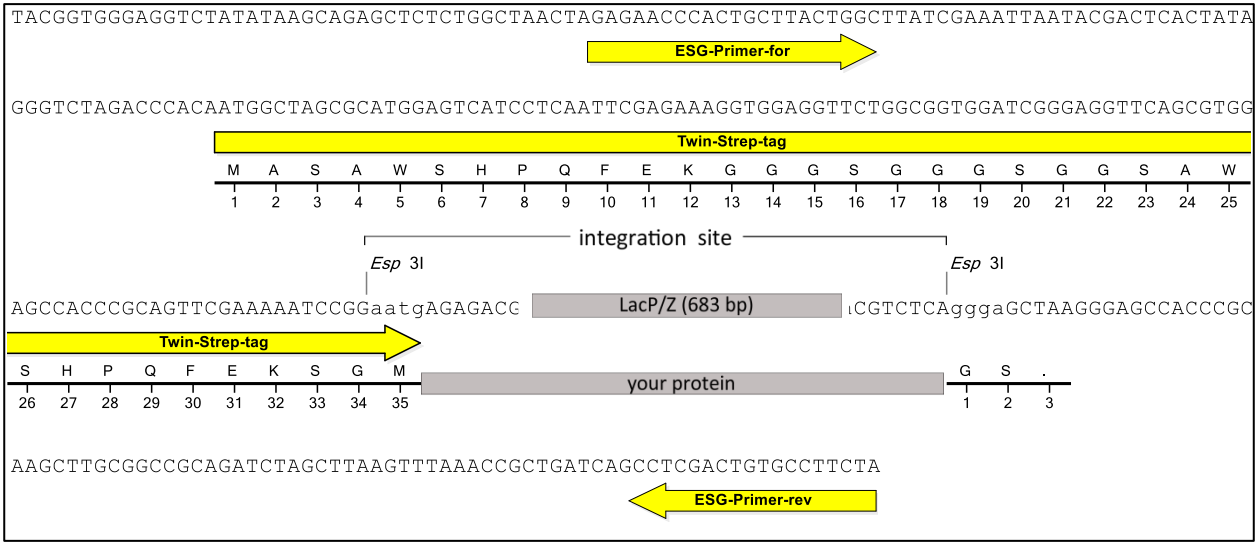

LacP/Z cassette = contains LacZ alpha fragment under control of a separate promoter, which allows alpha complementation of LacZ mutations such as LacZΔM15 as in *E. coli* DH5α or TOP10.

your protein = after StarGate cloning using *Esp3I* your gene of interest will be located here

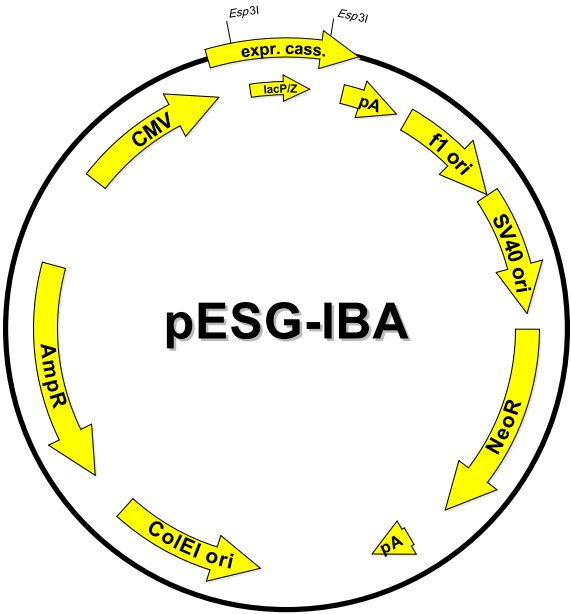

Features	from bp	to bp	Sequencing primer
f1 origin	259	687	ESG-Primer-for 5' – GAGAACCCACTGCTTACTGGC –3' ESG-Primer-rev 5' – TAGAAGGCACAGTCGAGG –3'
SV40 ori	692	1035	
Neomycin resistance gene	1097	1891	
ColEI ori	2637	3222	
Ampicillin resistance gene	3393	4253	
CMV promoter	4621	5208	
forward primer binding site	5221	5241	
Twin-Strep-tag®	5284	5388	
LacZ alpha fragment	5614	6015	
reverse primer binding site	6146	6163	
polyA signal sequence	1	213	
total vector length		6163	