



# SZABO SCANDIC

Part of Europa Biosite

## Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!  
See the following pages for more information!



### Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

### Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

### SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

[mail@szabo-scandic.com](mailto:mail@szabo-scandic.com)

[www.szabo-scandic.com](http://www.szabo-scandic.com)

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 

## FcR blocking Reagent

Reference	Size
FCR-2ML	200 test
FCR-4ML	400 test

### PRODUCT DESCRIPTION

FcR Blocking Reagent for human has high affinity receptor for IgG (Fc RI, CD64) but it is not an anti-CD64, which allows its use with pure antibodies.

Cells from patients treated with IFN-gamma or G-CSF or from patients subject to acute inflammation upregulate immunoglobulin Fc receptors (FcR) and are therefore particularly sensitive to this nonspecific Cy-binding effect.

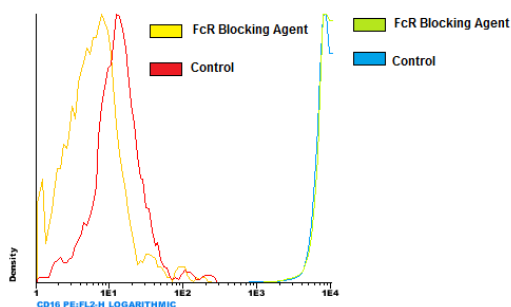


Fig. 1. Fluorescence profiles of normal human peripheral blood CD16 PE stained with FcR Blocking Reagent (yellow and green histogram) or without FcR Blocking Reagent (red and blue histogram).

**Recommended usage:** This product has been designed to block receptors constant fraction of monocytes and macrophages to avoid nonspecific binding of antibodies and fluorochromes and increasing the affinity of the antibodies to their targets which can complicate interpretation of flow cytometric results.

**Presentation:** liquid

**Storage instructions:** Store in the dark, refrigerated between 2 °C and 8 °C.

**Reagent provided:** FcR Blocking Reagent for human (10 µl/test) supplied in buffer containing protein stabilizer solution, PBS 20 mM and sodium azide 0,09%, pH 7,2.

**Recommendations and warnings:** The reagents contain sodium azide. In acid conditions, it is transformed into hydrazoic acid, a highly toxic compound. Avoid microbial contamination of the reagent. Protect from light.

Deviations from the recommended procedure could invalidate the analysis results.

For professional use only.

Reagents should not be used if any evidence of deterioration is observed. For more information, please contact our technical service: [tech@immunostep.com](mailto:tech@immunostep.com)

### FcR Blocking Reagent protocol for indirect labelling of human cells

1. Transfer the sample of blood to a 12 x 75 mm polystyrene test tube ( $10^5$  to  $10^8$  cells/test).
2. Add 10 µL of FcR Blocking Reagent (the optimal volume should be determined by the individual laboratory).
3. Mix well and incubate for 10 minutes at room temperature (20-25 °C).
4. Add purified antibodies according to manufacturer's recommendation and mix gently with a vortex mixer.
5. The recommended negative control is a non-reactive PU- antibody of the same isotype.
6. Incubate in the dark at room temperature at 4 °C for 30 minutes or at room temperature (20-25 °C) for 15 minutes.
7. Add 2 mL 0.01 mol/L PBS (It better that it containing 2% bovine serum albumin) and resuspend the cells by using a vortex mixer.
8. A washing is made with centrifugation at 2000 rpm for 5 min in order to remove the McAb not bound to its antigen.
9. A secondary conjugated with some fluorochrome is added and the mixture is incubated at room temperature for 15 min in the darkness. The absence of light is necessary so that the fluorochrome will not deteriorate since it shows a high degree of photoinstability.

10. After the incubation period, an erythrocyte-lysing solution is added at the amount recommended by the manufacturer and the mixture is incubated at room temperature in the darkness (the blood should be well mixed with the lysing solution).
  11. The tubes are centrifuged at 2000 rpm for 5 minutes. The supernatant is removed with a Pasteur pipette or with a vacuum pump.
  12. The cell pellet is resuspended and a final wash is made with 3-5 mL of PBS at 2000 rpm for 5 min.
  13. After removing the supernatant and resuspending the cell pellet, some 300  $\mu$ L of PBS is added and the readings on the flow cytometer are recorded.
  14. Analyse on a flow cytometer or store at 2-8 °C in the dark until analysis. Samples can be run up to 24 hours after lysis.
8. Centrifuge at 1000 x g for 5 minutes. Gently aspirate the supernatant and discard it leaving approximately 50  $\mu$ L of fluid.
  9. Add 2 mL 0.01 mol/L PBS (It better that it containing 2% bovine serum albumin) and resuspend the cells by using a vortex mixer.
  10. Centrifuge at 1000 x g for 5 minutes. Gently aspirate the supernatant and discard it leaving approximately 50  $\mu$ L of fluid.
  11. Resuspend pellet in an appropriate fluid for flow cytometry, e.g. 0.3 mL PBS. The PBS should contain 1% paraformaldehyde (fixative) if samples are not analysed the same day.
  12. Analyse on a flow cytometer or store at 2-8 °C in the dark until analysis. Samples can be run up to 24 hours after lysis.

#### **FcR Blocking Reagent protocol for direct labelling of human cells**

1. Transfer 100  $\mu$ L of anticoagulated (EDTA) blood to a 12 x 75 mm polystyrene test tube ( $10^6$  cells).
2. Add 10  $\mu$ L of FcR Blocking Reagent (the optimal volume should be determined by the individual laboratory).
3. Mix well and incubate for 10 minutes at room temperature (20-25 °C).
4. Add antibodies according to manufacturer's recommendation and mix gently with a vortex mixer.
5. The recommended negative control is a non-reactive conjugated antibody of the same isotype.
6. Incubate in the dark at room temperature at 4 °C for 30 minutes or at room temperature (20-25 °C) for 15 minutes.
7. Add 1,5 mL of Lysing Solution to each sample and mix gently with a vortex mixer. Incubate for 10 minutes at room temperature in the dark.

#### **WARRANTY**

Warranted only to conform to the quantity and contents stated on the label or in the product labelling at the time of delivery to the customer. Immunostep disclaims hereby other warranties. Immunostep's sole liability is limited to either the replacement of the products or refund of the purchase price.

#### **MANUFACTURED BY**



**Immunostep S.L.**  
Avda. Universidad de Coimbra, s/n  
Cancer Research Center (CIC)  
Campus Miguel de Unamuno  
37007 Salamanca (Spain)  
Tel. (+34) 923 294 827  
[www.immunostep.com](http://www.immunostep.com)