



SZABO SCANDIC

Part of Europa Biosite

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!
See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7


mail@szabo-scandic.com

www.szabo-scandic.com

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 

Dabsyl chloride

Cat. No.:	HY-101890
CAS No.:	56512-49-3
Molecular Formula:	C ₁₄ H ₁₄ ClN ₃ O ₂ S
Molecular Weight:	323.8
Target:	Fluorescent Dye
Pathway:	Others
Storage:	4°C, protect from light, stored under nitrogen * In solvent : -80°C, 6 months; -20°C, 1 month (protect from light, stored under nitrogen)

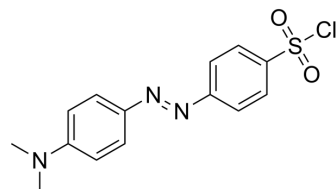

SOLVENT & SOLUBILITY

In Vitro	DMSO : 2 mg/mL (6.18 mM; Need ultrasonic)				
	Preparing Stock Solutions	Solvent Concentration Mass	1 mg	5 mg	10 mg
		1 mM	3.0883 mL	15.4416 mL	30.8833 mL
		5 mM	0.6177 mL	3.0883 mL	6.1767 mL
		10 mM	---	---	---
Please refer to the solubility information to select the appropriate solvent.					
In Vivo	1. Add each solvent one by one: 10% DMSO >> 40% PEG300 >> 5% Tween-80 >> 45% saline Solubility: ≥ 0.52 mg/mL (1.61 mM); Clear solution				
	2. Add each solvent one by one: 10% DMSO >> 90% (20% SBE-β-CD in saline) Solubility: ≥ 0.52 mg/mL (1.61 mM); Clear solution				

BIOLOGICAL ACTIVITY

Description	Dabsyl chloride is an amine derivatizing agent, able to give rise to stable products that can be easily monitored spectrophotometrically at 460 nm; Dabsyl chloride also used for labeling amino acids.
In Vitro	Dabsyl chloride can give rise to mono-Dabsyl and bis-Dabsyl derivatives in the presence of multiple amino groups. Furthermore with respect to OPA derivatization, Dabsyl chloride can react with primary and also with secondary amines ^[1] . MCE has not independently confirmed the accuracy of these methods. They are for reference only.

PROTOCOL

Cell Assay ^[1]

Selected mouse brain samples from either cortical or striatal regions (100 mg wet weight) and neuroblastoma cells (SH-SY5Y) pellet derived from 25 cm² flask are treated with 500 µL of 0.1 M HCl containing 0.2% TDGA, sonicated for 10 min (only for brain tissue), and then centrifuged at 14000 g for 30 min. The supernatant is freeze-dried. 50 µL of reaction buffer and 100 µL of 15 mM Dabsyl chloride are added to the tube and derivatized.

MCE has not independently confirmed the accuracy of these methods. They are for reference only.

CUSTOMER VALIDATION

- Lwt-Food Sci Technol. 15 January 2022, 112783.
- Int J Food Sci Tech. 2020 Aug 11.

See more customer validations on www.MedChemExpress.com

REFERENCES

[1]. Francioso A, et al. HPLC Determination of Bioactive Sulfur Compounds, Amino Acids and Biogenic Amines in Biological Specimens. Adv Exp Med Biol. 2017;975:535-549.

Caution: Product has not been fully validated for medical applications. For research use only.

Tel: 609-228-6898

Fax: 609-228-5909

E-mail: tech@MedChemExpress.com

Address: 1 Deer Park Dr, Suite Q, Monmouth Junction, NJ 08852, USA