



SZABO SCANDIC

Part of Europa Biosite

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!
See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7


mail@szabo-scandic.com

www.szabo-scandic.com

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 

VU6005649

Cat. No.:	HY-107982
CAS No.:	2137047-43-7
Molecular Formula:	C ₁₆ H ₁₂ F ₅ N ₃ O
Molecular Weight:	357.28
Target:	mGluR
Pathway:	GPCR/G Protein; Neuronal Signaling
Storage:	<div> <div>Powder</div> <div>-20°C 3 years</div> <div>4°C 2 years</div> </div> <div> <div>In solvent</div> <div>-80°C 2 years</div> <div>-20°C 1 year</div> </div>

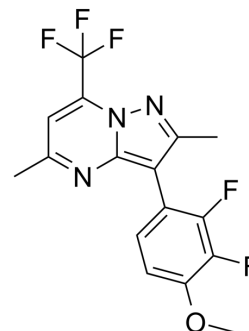

SOLVENT & SOLUBILITY

In Vitro	DMSO : 50 mg/mL (139.95 mM; Need ultrasonic)					
	Preparing Stock Solutions	<div><div>Solvent</div><div>Concentration</div></div>	Mass	1 mg	5 mg	10 mg
		1 mM		2.7989 mL	13.9946 mL	27.9893 mL
		5 mM		0.5598 mL	2.7989 mL	5.5979 mL
		10 mM		0.2799 mL	1.3995 mL	2.7989 mL
Please refer to the solubility information to select the appropriate solvent.						
In Vivo	1. Add each solvent one by one: 10% DMSO >> 40% PEG300 >> 5% Tween-80 >> 45% saline Solubility: ≥ 2.5 mg/mL (7.00 mM); Clear solution					
	2. Add each solvent one by one: 10% DMSO >> 90% corn oil Solubility: ≥ 2.5 mg/mL (7.00 mM); Clear solution					

BIOLOGICAL ACTIVITY

Description	VU6005649 is a CNS penetrant mGlu _{7/8} receptor agonist with EC ₅₀ s of 0.65 μM and 2.6 μM for mGlu ₇ receptor and mGlu ₈ receptor, respectively.	
IC ₅₀ & Target	mGlu7 Receptor 0.65 μM (EC ₅₀)	mGlu8 Receptor 2.6 μM (EC ₅₀)
In Vitro	<p>VU6005649 is a CNS penetrant mGlu_{7/8} receptor agonist with EC₅₀s of 0.65 μM and 2.6 μM for mGlu₇ receptor and mGlu₈ receptor, respectively. VU6005649 displays a terminal K_p of 2.43 with total brain levels ~9× above the mGlu₇ positive allosteric modulator (PAM) in vitro EC₅₀^[1].</p> <p>MCE has not independently confirmed the accuracy of these methods. They are for reference only.</p>	

In Vivo

When VU6005649 (compound 9f) is dosed at 30 mg/kg IP in 10% Tween 80/H₂O (0.75 mg/kg. s.c. amphetamine), no efficacy is observed in this assay. VU6005649 shows modest but significant pro-cognitive effects on associative learning in wild-type mice and the first example of efficacy of an mGlu_{7/8} positive allosteric modulator (PAM) in this model^[1]. MCE has not independently confirmed the accuracy of these methods. They are for reference only.

PROTOCOL

Animal Administration ^[1]

Tissue distribution studies with VU6005649 (compound 9f) in mice are performed by formulating VU6005649 in 10% polysorbate 80 and dosing via intraperitoneal injection to 20 week old female C57/Bl6 mice (3 per time point). At 0.25, 0.5, 1, 3, and 6 hours post dose, animals are euthanized and decapitated, blood is collected via cardiac puncture and the brains are removed, thoroughly washed in cold phosphate-buffered saline, and immediately frozen on dry ice^[1]. MCE has not independently confirmed the accuracy of these methods. They are for reference only.

REFERENCES

[1]. Abe M, et al. Discovery of VU6005649, a CNS Penetrant mGlu_{7/8} Receptor PAM Derived from a Series of Pyrazolo[1,5-a]pyrimidines. ACS Med Chem Lett. 2017 Sep 1;8(10):1110-1115.

Caution: Product has not been fully validated for medical applications. For research use only.

Tel: 609-228-6898

Fax: 609-228-5909

E-mail: tech@MedChemExpress.com

Address: 1 Deer Park Dr, Suite Q, Monmouth Junction, NJ 08852, USA