



SZABO SCANDIC

Part of Europa Biosite

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!
See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7


mail@szabo-scandic.com

www.szabo-scandic.com

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 

ETP-45658

Cat. No.:	HY-110109
CAS No.:	1198357-79-7
Molecular Formula:	C ₁₆ H ₁₇ N ₅ O ₂
Molecular Weight:	311.34
Target:	PI3K; DNA-PK; mTOR
Pathway:	PI3K/Akt/mTOR; Cell Cycle/DNA Damage
Storage:	Powder -20°C 3 years 4°C 2 years In solvent -80°C 6 months -20°C 1 month

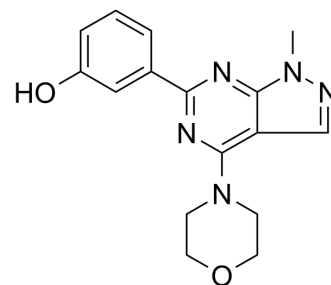

SOLVENT & SOLUBILITY

In Vitro	DMSO : 250 mg/mL (802.98 mM; Need ultrasonic)					
	Preparing Stock Solutions	Solvent Concentration	Mass	1 mg	5 mg	10 mg
		1 mM		3.2119 mL	16.0596 mL	32.1192 mL
		5 mM		0.6424 mL	3.2119 mL	6.4238 mL
		10 mM		0.3212 mL	1.6060 mL	3.2119 mL
Please refer to the solubility information to select the appropriate solvent.						
In Vivo	1. Add each solvent one by one: 10% DMSO >> 40% PEG300 >> 5% Tween-80 >> 45% saline Solubility: ≥ 2.08 mg/mL (6.68 mM); Clear solution					
	2. Add each solvent one by one: 10% DMSO >> 90% (20% SBE-β-CD in saline) Solubility: ≥ 2.08 mg/mL (6.68 mM); Clear solution					
	3. Add each solvent one by one: 10% DMSO >> 90% corn oil Solubility: ≥ 2.08 mg/mL (6.68 mM); Clear solution					

BIOLOGICAL ACTIVITY

Description	ETP-45658 is a potent PI3K inhibitor, with IC ₅₀ s of 22.0 nM, 39.8 nM, 129.0 nM and 717.3 nM for PI3Kα, PI3Kδ, PI3Kβ and PI3Kγ, respectively. ETP-45658 also can inhibit DNA-PK (IC ₅₀ =70.6 nM) and mTOR (IC ₅₀ =152.0 nM). ETP-45658 can be used for the research of cancer ^{[1][2]} .			
IC ₅₀ & Target	PI3Kα 22.0 nM (IC ₅₀)	PI3Kδ 39.8 nM (IC ₅₀)	PI3Kβ 129.0 nM (IC ₅₀)	PI3Kγ 717.3 nM (IC ₅₀)
	DNA-PK	mTOR		

	70.6 nM (IC ₅₀)	152.0 nM (IC ₅₀)
In Vitro	ETP-45658 (10 μM; 4 h) decreases in the phosphorylation of FOXO3a, Gsk3-β and p70 S6K in U2OS cells^[1]. ETP-45658 inhibits the proliferation of MCF7, PC3, 786-O, HTC116, and U251 cells, with EC₅₀s of 0.48 μM, 0.49 μM, 2.62 μM, 3.53 μM, and 5.56 μM, respectively^[1]. ETP-45658 (10 μM; 24 h) induces a clear G1 arrest of PC3 cells^[1]. ETP-45658 inhibits the mutant PI3Kα proteins, H1047R and E545K, with IC₅₀s of 16.8 nM and 13.1 nM, respectively^[1]. ETP-45658 (5 nM-11.1 μM; 1 h) induces a dose-dependent increase of GFP-FOXO nuclear translocation in U2foxRELOC cells^[1]. ETP-45658 (10 μM; 1 h) decrease the expression of cyclin D1 and p-Akt on serine 473 in U2OS cells^[1]. MCE has not independently confirmed the accuracy of these methods. They are for reference only. Western Blot Analysis^[1]	
	Cell Line:	U2OS cells
	Concentration:	10 μM
	Incubation Time:	4 hours
	Result:	Caused a 95% decrease in the phosphorylation of FOXO3a at threonine 32 and a 55% reduction on Gsk3-β at serine 9.
		Reduced the phosphorylation of the mTOR substrate p70 S6K.
In Vivo	ETP-45658 (22.7 mg/kg) decreases the level of phosphorylated Akt on serine 473 in the mammary ducts of transgenic mice^[1]. MCE has not independently confirmed the accuracy of these methods. They are for reference only.	

REFERENCES

- [1]. Link W, et, al. Chemical interrogation of FOXO3a nuclear translocation identifies potent and selective inhibitors of phosphoinositide 3-kinases. J Biol Chem. 2009 Oct 9;284(41):28392-28400.
- [2]. Hill R, et, al. A novel phosphatidylinositol 3-kinase (PI3K) inhibitor directs a potent FOXO-dependent, p53-independent cell cycle arrest phenotype characterized by the differential induction of a subset of FOXO-regulated genes. Breast Cancer Res. 2014 Dec 9;16(6):482.

Caution: Product has not been fully validated for medical applications. For research use only.

Tel: 609-228-6898

Fax: 609-228-5909

E-mail: tech@MedChemExpress.com

Address: 1 Deer Park Dr, Suite Q, Monmouth Junction, NJ 08852, USA