



SZABO SCANDIC

Part of Europa Biosite

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!
See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7


mail@szabo-scandic.com

www.szabo-scandic.com

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 

Ethylidene-4-nitrophenyl- α -D-Maltoheptaoside

Cat. No.:	HY-112835
CAS No.:	96597-16-9
Molecular Formula:	C ₅₀ H ₇₇ NO ₃₈
Molecular Weight:	1300.13
Target:	Amylases
Pathway:	Metabolic Enzyme/Protease
Storage:	4°C, protect from light, stored under nitrogen * In solvent : -80°C, 6 months; -20°C, 1 month (protect from light, stored under nitrogen)

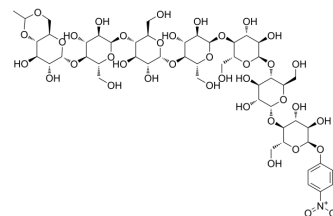

SOLVENT & SOLUBILITY

In Vitro	DMSO : 100 mg/mL (76.92 mM; Need ultrasonic)				
	Preparing Stock Solutions	Solvent Concentration Mass	1 mg	5 mg	10 mg
		1 mM	0.7692 mL	3.8458 mL	7.6915 mL
		5 mM	0.1538 mL	0.7692 mL	1.5383 mL
		10 mM	0.0769 mL	0.3846 mL	0.7692 mL
Please refer to the solubility information to select the appropriate solvent.					
In Vivo	1. Add each solvent one by one: 10% DMSO >> 40% PEG300 >> 5% Tween-80 >> 45% saline Solubility: ≥ 2.5 mg/mL (1.92 mM); Clear solution				
	2. Add each solvent one by one: 10% DMSO >> 90% (20% SBE-β-CD in saline) Solubility: ≥ 2.5 mg/mL (1.92 mM); Clear solution				

BIOLOGICAL ACTIVITY

Description	Ethylidene-4-nitrophenyl- α -D-Maltoheptaoside (pNP-G7) is the substrate of α -amylase. Ethylidene-4-nitrophenyl- α -D-Maltoheptaoside can be used to measure α -amylase activity ^[1] .
-------------	--

REFERENCES

- [1]. Lorentz K. Routine alpha-amylase assay using protected 4-nitrophenyl-1, 4-alpha-D-maltoheptaoside and a novel alpha-glucosidase. Clin Chem. 2000 May;46(5):644-9.
- [2]. Lorentz K. Routine alpha-amylase assay using protected 4-nitrophenyl-1, 4-alpha-D-maltoheptaoside and a novel alpha-glucosidase. Clin Chem. 2000 May;46(5):644-9.

Caution: Product has not been fully validated for medical applications. For research use only.

Tel: 609-228-6898

Fax: 609-228-5909

E-mail: tech@MedChemExpress.com

Address: 1 Deer Park Dr, Suite Q, Monmouth Junction, NJ 08852, USA