



SZABO SCANDIC

Part of Europa Biosite

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!
See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC Handels GmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7


mail@szabo-scandic.com

www.szabo-scandic.com

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 

GSK126

Cat. No.:	HY-13470
CAS No.:	1346574-57-9
Molecular Formula:	C ₃₁ H ₃₈ N ₆ O ₂
Molecular Weight:	526.67
Target:	Histone Methyltransferase
Pathway:	Epigenetics
Storage:	<div>Powder</div> <div>-20°C 3 years</div> <div>4°C 2 years</div> <div>In solvent</div> <div>-80°C 2 years</div> <div>-20°C 1 year</div>

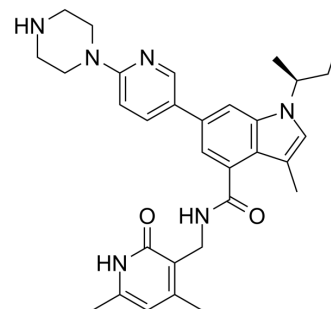

SOLVENT & SOLUBILITY

In Vitro

DMSO : 12.5 mg/mL (23.73 mM; Need ultrasonic)

	Solvent Concentration	Mass	1 mg	5 mg	10 mg
Preparing Stock Solutions	1 mM		1.8987 mL	9.4936 mL	18.9872 mL
	5 mM		0.3797 mL	1.8987 mL	3.7974 mL
	10 mM		0.1899 mL	0.9494 mL	1.8987 mL

Please refer to the solubility information to select the appropriate solvent.

In Vivo

- Add each solvent one by one: 20% SBE-β-CD adjusted to pH 4-4.5 with 1 N acetic
Solubility: 20 mg/mL (37.97 mM); Clear solution; Need ultrasonic and warming
- Add each solvent one by one: 50% PEG300 >> 50% saline
Solubility: 10 mg/mL (18.99 mM); Suspended solution; Need ultrasonic
- Add each solvent one by one: 10% DMSO >> 40% PEG300 >> 5% Tween-80 >> 45% saline
Solubility: ≥ 1.25 mg/mL (2.37 mM); Clear solution
- Add each solvent one by one: 10% DMSO >> 90% (20% SBE-β-CD in saline)
Solubility: 1.25 mg/mL (2.37 mM); Suspended solution; Need ultrasonic
- Add each solvent one by one: 10% DMSO >> 90% corn oil
Solubility: ≥ 1.25 mg/mL (2.37 mM); Clear solution
- Add each solvent one by one: 5% DMSO >> 40% PEG300 >> 5% Tween-80 >> 50% saline
Solubility: ≥ 0.82 mg/mL (1.56 mM); Clear solution
- Add each solvent one by one: 5% DMSO >> 95% (20% SBE-β-CD in saline)
Solubility: ≥ 0.82 mg/mL (1.56 mM); Clear solution

BIOLOGICAL ACTIVITY

Description	GSK126 (GSK2816126A) is a potent, highly selective inhibitor of EZH2 methyltransferase with an IC ₅₀ of 9.9 nM ^[1] .
IC₅₀ & Target	EZH2 9.9 nM (IC ₅₀)
In Vitro	<p>GSK126 potently inhibits both wild-type and mutant EZH2 methyltransferase activity with similar potencies (K_i=0.5-3 nM) independent of substrate used, and is competitive with S-adenosyl-methionine (SAM) and non-competitive with peptide substrates. GSK126 is highly selective against other methyltransferases and multiple other protein classes (EZH1, IC₅₀=680 nM)^[1]. Treatment of three SCLC cell lines with GSK126, induces growth inhibition. SCLC cell lines (Lu130, H209, and DMS53) are treated with 0.5, 2, and 8 μM GSK126, and growth curve is analyzed by WST-8 assay. Inhibition of cellular growth by GSK126 treatment is observed at 8 μM in all the three cell lines, while Lu130 and H209 are more sensitive to GSK126, even at lower doses^[2].</p> <p>MCE has not independently confirmed the accuracy of these methods. They are for reference only.</p>
In Vivo	<p>GSK126 is administered intraperitoneally at a dose volume of 0.2 mL per 20 g body weight in female beige SCID mice. GSK126 effectively inhibits the proliferation of EZH2 mutant DLBCL cell lines and markedly inhibits the growth of EZH2 mutant DLBCL xenografts in mice^[1].</p> <p>MCE has not independently confirmed the accuracy of these methods. They are for reference only.</p>

PROTOCOL

Kinase Assay ^[1]	<p>The five-member PRC2 complex (Flag-EZH2, EED, SUZ12, AEBP2, RbAp48) containing either wild-type or mutant (A677G, Y641N, Y641C, Y641H, Y641S or Y641F) EZH2 is prepared. GSK126 is dissolved in DMSO and tested at concentrations of 0.6 nM to 300 nM with a final DMSO concentration of 2.5%. In contrast to wild-type EZH2 which prefers H3K27me0 as a substrate in vitro, EZH2 Y641 mutants prefer H3K27me2 and have little activity with H3K27me0 or H3K27me1. The A677G mutant is distinct from both the wild-type and Y641 mutant forms of EZH2 in that it efficiently methylates H3K27me0, H3K27me1, and H3K27me2; therefore, histone H3 peptides (residues 21-44; 10 μM final) with either K27me0 (wild type, A677G EZH2), K27me1 (A677G EZH2), or K27me2 (A677G, Y641N, Y641C, Y641H, Y641S and Y641F EZH2) are used as methyltransferase substrates. GSK126 is added to plates followed by addition of 6 nM EZH2 complex and peptide. As the potency of GSK126 is at or near the tight binding limit of an assay run at [SAM]=K_m, IC₅₀ values are measured at a high concentration of the competitive substrate SAM relative to its K_m (7.5 μM SAM where the SAM K_m is 0.3 μM)^[1].</p> <p>MCE has not independently confirmed the accuracy of these methods. They are for reference only.</p>
Cell Assay ^[2]	<p>JUB- and PTRF-introduced DMS53 cells are seeded at density of 1×10³ cells/well in 96-well plate, and cellular growth is analyzed using WST-8 kit at 12, 36, 60, and 84 h. Cellular growth of Lu130, H209, and DMS53 with treatment by DZNep or GSK126 is also analyzed using WST-8 kit. DZNep is dissolved in PBS at 5 mM, and cells are cultured at the final concentration of 5 μM. GSK126 is dissolved in DMSO at 10 mM, and cells are cultured at 0.5, 2, and 8 μM^[2].</p> <p>MCE has not independently confirmed the accuracy of these methods. They are for reference only.</p>
Animal Administration ^[1]	<p>Mice^[1]</p> <p>GSK126 or vehicle is administered intraperitoneally at a dose volume of 0.2 mL per 20 g body weight. Pfeiffer or KARPAS-422 cells (1×10⁷) in 100% Matrigel are implanted subcutaneously in female beige SCID mice. Tumors are measured with calipers, and block randomized according to tumour size into treatment groups. For efficacy studies, 10 mice are randomized in each treatment group before the initiation of dosing and GSK126 treatment is initiated once the tumour volumes are approximately 200 mm³ in the Pfeiffer and KARPAS-422 studies and 500 mm³ in the KARPAS-422 intermittent dosing study. Mice are weighed and tumors measured with calipers twice weekly. Two-tailed t-tests are conducted assuming two samples of equal variance.</p> <p>MCE has not independently confirmed the accuracy of these methods. They are for reference only.</p>

CUSTOMER VALIDATION

-
- Cell. 2018 Sep 20;175(1):186-199.e19.
 - Blood. 2022 Aug 19;blood.2022015668.
 - Mol Cell. 2023 Nov 20:S1097-2765(23)00914-0.
 - J Exp Med. 2018 May 7;215(5):1365-1382.
 - Cell Death Differ. 2022 May 14.

See more customer validations on www.MedChemExpress.com

REFERENCES

- [1]. McCabe MT, et al. EZH2 inhibition as a therapeutic strategy for lymphoma with EZH2-activating mutations. Nature. 2012 Dec 6;492(7427):108-12.
- [2]. Sato T, et al. PRC2 overexpression and PRC2-target gene repression relating to poorer prognosis in small cell lung cancer. Sci Rep. 2013 May 29;3:1911.
-

Caution: Product has not been fully validated for medical applications. For research use only.

Tel: 609-228-6898

Fax: 609-228-5909

E-mail: tech@MedChemExpress.com

Address: 1 Deer Park Dr, Suite Q, Monmouth Junction, NJ 08852, USA