



SZABO SCANDIC

Part of Europa Biosite

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!
See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7


mail@szabo-scandic.com

www.szabo-scandic.com

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 

AG-1478

Cat. No.:	HY-13524		
CAS No.:	153436-53-4		
Molecular Formula:	C ₁₆ H ₁₄ ClN ₃ O ₂		
Molecular Weight:	316		
Target:	EGFR; Influenza Virus; HCV		
Pathway:	JAK/STAT Signaling; Protein Tyrosine Kinase/RTK; Anti-infection		
Storage:	Powder	-20°C	3 years
		4°C	2 years
	In solvent	-80°C	2 years
		-20°C	1 year

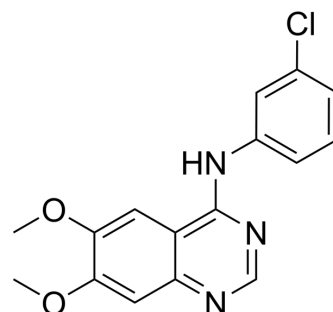

SOLVENT & SOLUBILITY

In Vitro	DMSO : 50 mg/mL (158.23 mM; Need ultrasonic)					
	Preparing Stock Solutions	<div><div>Solvent</div><div>Concentration</div></div>	Mass	1 mg	5 mg	10 mg
		1 mM		3.1646 mL	15.8228 mL	31.6456 mL
		5 mM		0.6329 mL	3.1646 mL	6.3291 mL
		10 mM		0.3165 mL	1.5823 mL	3.1646 mL
Please refer to the solubility information to select the appropriate solvent.						
In Vivo	1. Add each solvent one by one: 1% CMC >> 0.5% Tween-80 Solubility: 5 mg/mL (15.82 mM); Suspended solution; Need ultrasonic					
	2. Add each solvent one by one: 10% DMSO >> 40% PEG300 >> 5% Tween-80 >> 45% saline Solubility: ≥ 2.08 mg/mL (6.58 mM); Clear solution					
	3. Add each solvent one by one: 10% DMSO >> 90% (20% SBE-β-CD in saline) Solubility: ≥ 2.08 mg/mL (6.58 mM); Clear solution					
	4. Add each solvent one by one: 10% DMSO >> 90% corn oil Solubility: ≥ 2.08 mg/mL (6.58 mM); Clear solution					

BIOLOGICAL ACTIVITY

Description	AG-1478 (Tyrphostin AG-1478) is a selective EGFR tyrosine kinase inhibitor with IC ₅₀ of 3 nM. AG-1478 has antiviral effects against HCV and encephalomyocarditis virus (EMCV).		
IC ₅₀ & Target	EGFR 3 nM (IC ₅₀)	HCV	EMCV

In Vitro	<p>AG-1478 (AG1478) is irreversible for growth regulation of human lung (A549) and prostate (DU145) cancer cell lines, cultured in chemically defined DMEM/F12 medium. AG-1478 seems to be more effective at lower concentrations, but is unable to completely inhibit growth of A549 cells^[1]. Inhibition of EGFR by specific tyrosine kinase inhibitor AG-1478 (AG1478) significantly decreases the angiotensin II-mediated synthesis of TGF-β and fibronectin by cardiac fibroblasts. EGFR is pharmacologically inhibited by small-molecule inhibitor AG-1478 with IC₅₀ of 4 nM^[2]. Both Polyfect (PF) and Superfect (SF) treatment lead to increased apoptosis in HEK 293 cells to a similar extent as assessed by flow cytometry. The antioxidant, tempol, significantly reduced dendrimer-mediated apoptosis for both PF and SF. AG-1478 (AG1478), at a 10-fold higher dose (100 μM) than used in signaling studies, is used as a positive control and significantly induced apoptosis in HEK 293 cells^[3]. MCE has not independently confirmed the accuracy of these methods. They are for reference only.</p>
In Vivo	<p>Administration of AG-1478 (AG1478) significantly reduces myocardial inflammation, fibrosis, apoptosis, and dysfunction in both two obese mouse models. ApoE^{-/-} mice are first fed with HFD for 8 weeks (ApoE-HFD), and then administrated with AG-1478 (10 mg/kg/day) or 542 (10 mg/kg/day) for another 8 weeks by oral gavage. AG-1478 or 542 treatment blocks HFD induced cardiac EGFR phosphorylation in vivo, without affecting the plasma level of low density lipoprotein (LDL) and total triglyceride (TG)^[2]. Administration of EGF (10 nM) leads to a robust and reproducible elevation in EGFR phosphorylation that can be blocked by AG-1478 (AG1478), a known inhibitor of EGFR phosphorylation. Increasing doses of Polyfect (PF) result in a significant reduction in EGF-induced EGFR phosphorylation (p<0.05) but this is to a lesser extent than observed with AG1478^[3]. MCE has not independently confirmed the accuracy of these methods. They are for reference only.</p>

PROTOCOL

Cell Assay ^[1]	<p>DU145 (HTB-81) and A549 (CCL-185) cells are seeded on 96-well plates at concentrations of 4×10³ cells/well in MEM (DU145 cells) or DMEM (A549 cells). Following 24 h of incubation, the culture medium is replaced with serum-free DMEM/F12 (1:1) supplemented with Transferrin (5 mg/mL), sodium selenite (2 ng/mL) and albumin (0.5 mg/mL) [DMEM/F12+]. After an additional 24 h of incubation (Day 0), the medium is replaced by serum-free DMEM/F12+ medium containing tyrosine kinase inhibitors: AG494, AG-1478 respectively in concentration ranges 1-20 μM and 0.1-8 μM. The incubation is continued for the next 24 h at 37°C in a humidified atmosphere. The modified crystal violet staining method (CV) and MTT assay are used to determine the influence of the tyrphostins on the proliferation of target cells. The absorbance is measured using a Tecan multiscan plate recorder^[1]. MCE has not independently confirmed the accuracy of these methods. They are for reference only.</p>
Animal Administration ^{[2][3]}	<p>Mice^[2] 28 C57BL/6 or ApoE^{-/-} mice are randomly divided into four weight-matched groups. 7 mice are fed with standard animal low-fat diet containing 10 kcal.% fat, 20 kcal.% protein and 70 kcal.% carbohydrate serve as a normal control group (Control or ApoE-LF), while the remaining 21 mice are fed with high-fat diet containing 60 kcal.% fat, 20 kcal.% protein and 20 kcal.% carbohydrate for 16 weeks. Since 9th week, AG-1478 or 542 are given daily by oral gavage at a dose of 10 mg/kg/day for the next 8 weeks. Mice in the Control and HFD groups are gavaged with vehicle (1% CMC-Na solution) only. At the day before the sacrifice of ApoE^{-/-} mice, doppler analysis is performed to determine the pathologic cardiac hypertrophy.</p> <p>Rats^[3] Male Wistar rats weighing about 300g are used in this study and divided into the following groups (N=5). Group 1: Non-diabetic (Control, C) animals, Group 2: C+PF (10mg/kg administered as a single intraperitoneal (i.p) injection) Group 3: C+SF (10mg/kg i.p); Group 4: C+AG-1478 (1 mg/kg i.p). Group 5: Rats bearing 4 weeks of diabetes (D) induced by a single i.p. injection of streptozotocin (55 mg/kg body weight); Group 6: D+PF (10 mg/kg i.p) Group 7: D+SF (10 mg/kg i.p) and Group 8: D+AG-1478 (1 mg/kg i.p). AG-1478 and dendrimer treatments are administered as single dose for 24h prior to sacrifice. Rat body weight and basal glucose levels are assessed before and after treatments just before sacrificing the animals. An automated blood glucose analyzer is used to assess blood glucose concentrations and rats with a blood glucose concentration above 250 mg/dL (approx. 14 mM) are declared diabetic as in previous studies. MCE has not independently confirmed the accuracy of these methods. They are for reference only.</p>

CUSTOMER VALIDATION

- Adv Sci (Weinh). 2023 Apr 14;e2204824.
- Cell Death Dis. 2021 May 18;12(6):509.
- Cell Death Dis. 2020 Jun 15;11(6):459.
- Cell Death Dis. 2019 Sep 9;10(9):649.
- Oncogene. 2021 Jul;40(30):4884-4893.

See more customer validations on www.MedChemExpress.com

REFERENCES

- [1]. Bojko A, et al. The effect of tyrphostins AG494 and AG1478 on the autocrine growth regulation of A549 and DU145 cells. *Folia Histochem Cytobiol.* 2012 Jul 5;50(2):186-95.
- [2]. Li W, et al. EGFR Inhibition Blocks Palmitic Acid-induced inflammation in cardiomyocytes and Prevents Hyperlipidemia-induced Cardiac Injury in Mice. *Sci Rep.* 2016 Apr 18;6:24580.
- [3]. Akhtar S, et al. Cationic Polyamidoamine Dendrimers as Modulators of EGFR Signaling In Vitro and In Vivo. *PLoS One.* 2015 Jul 13;10(7):e0132215.
- [4]. Dorobantu CM, et al. Tyrphostin AG1478 Inhibits Encephalomyocarditis Virus and Hepatitis C Virus by Targeting Phosphatidylinositol 4-Kinase III α . *Antimicrob Agents Chemother.* 2016 Sep 23;60(10):6402-6.

Caution: Product has not been fully validated for medical applications. For research use only.

Tel: 609-228-6898

Fax: 609-228-5909

E-mail: tech@MedChemExpress.com

Address: 1 Deer Park Dr, Suite Q, Monmouth Junction, NJ 08852, USA