



SZABO SCANDIC

Part of Europa Biosite

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!
See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7


mail@szabo-scandic.com

www.szabo-scandic.com

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 

Chromenone 1

Cat. No.: HY-143891
CAS No.: 1639929-29-5
Molecular Formula: C₁₈H₁₀F₃N₃O₂
Molecular Weight: 357.29
Target: TGF-β Receptor
Pathway: TGF-beta/Smad
Storage: 4°C, protect from light
 * In solvent : -80°C, 6 months; -20°C, 1 month (protect from light)

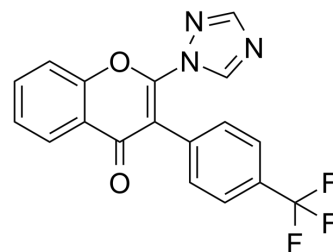

SOLVENT & SOLUBILITY

In Vitro	DMSO : 100 mg/mL (279.88 mM; Need ultrasonic)					
	Preparing Stock Solutions	Solvent Concentration	Mass	1 mg	5 mg	10 mg
		1 mM	2.7988 mL	13.9942 mL	27.9885 mL	
		5 mM	0.5598 mL	2.7988 mL	5.5977 mL	
		10 mM	0.2799 mL	1.3994 mL	2.7988 mL	
Please refer to the solubility information to select the appropriate solvent.						
In Vivo	1. Add each solvent one by one: 10% DMSO >> 40% PEG300 >> 5% Tween-80 >> 45% saline Solubility: ≥ 2.5 mg/mL (7.00 mM); Clear solution					
	2. Add each solvent one by one: 10% DMSO >> 90% corn oil Solubility: ≥ 2.5 mg/mL (7.00 mM); Clear solution					

BIOLOGICAL ACTIVITY

Description	Chromenone 1 is a potent osteogenic bone morphogenetic protein (BMP) potentiator. Chromenone 1 exhibits a unique mode of action as it induces a pronounced, kinase-independent, negative TGFβ feedback that enhances nuclear BMP-Smad signaling outputs ^[1] .
In Vitro	Chromenone 1 rescues DMH-1-induced cardiogenesis (IC ₅₀ =0.2 μM) and efficiently stimulates BMP-dependent osteogenesis from C2C12 myoblasts ^[1] . MCE has not independently confirmed the accuracy of these methods. They are for reference only.

REFERENCES

Caution: Product has not been fully validated for medical applications. For research use only.

Tel: 609-228-6898

Fax: 609-228-5909

E-mail: tech@MedChemExpress.com

Address: 1 Deer Park Dr, Suite Q, Monmouth Junction, NJ 08852, USA