



SZABO SCANDIC

Part of Europa Biosite

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!
See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7


mail@szabo-scandic.com

www.szabo-scandic.com

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 

DLin-M-K-DMA

Cat. No.:	HY-154803
CAS No.:	1221271-55-1
Molecular Formula:	C ₄₂ H ₇₇ NO ₂
Molecular Weight:	628.07
Target:	Others
Pathway:	Others
Storage:	Pure form -20°C 3 years 4°C 2 years In solvent -80°C 6 months -20°C 1 month

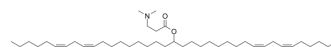

SOLVENT & SOLUBILITY

In Vitro	DMSO : 100 mg/mL (159.22 mM; Need ultrasonic)					
	Preparing Stock Solutions	Solvent Concentration	Mass	1 mg	5 mg	10 mg
		1 mM		1.5922 mL	7.9609 mL	15.9218 mL
		5 mM		0.3184 mL	1.5922 mL	3.1844 mL
		10 mM		0.1592 mL	0.7961 mL	1.5922 mL
Please refer to the solubility information to select the appropriate solvent.						
In Vivo	1. Add each solvent one by one: 10% DMSO >> 40% PEG300 >> 5% Tween-80 >> 45% saline Solubility: ≥ 2.5 mg/mL (3.98 mM); Clear solution					
	2. Add each solvent one by one: 10% DMSO >> 90% (20% SBE-β-CD in saline) Solubility: 2.5 mg/mL (3.98 mM); Suspended solution; Need ultrasonic					
	3. Add each solvent one by one: 10% DMSO >> 90% corn oil Solubility: ≥ 2.5 mg/mL (3.98 mM); Clear solution					

BIOLOGICAL ACTIVITY

Description	DLin-M-K-DMA (DLin-M-C2-DMA) is a cationic lipid that can be used in the construction of biological delivery vectors ^[1] .
-------------	---

REFERENCES

[1]. Michael J, et al. Improved amino lipids and methods for the delivery of nucleic acids. WO2010042877A1.

Caution: Product has not been fully validated for medical applications. For research use only.

Tel: 609-228-6898

Fax: 609-228-5909

E-mail: tech@MedChemExpress.com

Address: 1 Deer Park Dr, Suite Q, Monmouth Junction, NJ 08852, USA