



SZABO SCANDIC

Part of Europa Biosite

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!
See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7


mail@szabo-scandic.com

www.szabo-scandic.com

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 

Hoechst 33342 analog

Cat. No.:	HY-15627
CAS No.:	178481-68-0
Molecular Formula:	C ₃₂ H ₃₇ Cl ₂ N ₇
Molecular Weight:	590.59
Target:	Fluorescent Dye
Pathway:	Others
Storage:	-20°C, protect from light * In solvent : -80°C, 6 months; -20°C, 1 month (protect from light)

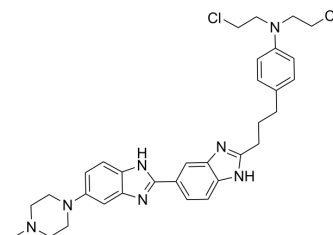

SOLVENT & SOLUBILITY

In Vitro	DMSO : 100 mg/mL (169.32 mM; Need ultrasonic)					
	Preparing Stock Solutions	Solvent Concentration	Mass	1 mg	5 mg	10 mg
		1 mM		1.6932 mL	8.4661 mL	16.9322 mL
		5 mM		0.3386 mL	1.6932 mL	3.3864 mL
		10 mM		0.1693 mL	0.8466 mL	1.6932 mL
Please refer to the solubility information to select the appropriate solvent.						
In Vivo	1. Add each solvent one by one: 10% DMSO >> 40% PEG300 >> 5% Tween-80 >> 45% saline Solubility: ≥ 2.5 mg/mL (4.23 mM); Clear solution					

BIOLOGICAL ACTIVITY

Description	Hoechst 33342 analog is a marker dye in Hoechst series. Hoechst is A live nuclear marker dye. Hoechst binds to the grooves in the DNA double strand, which tends to be A/ T-rich DNA strand. Although it binds to all nucleic acids, the A/ T-rich double strand DNA significantly enhances fluorescence intensity Therefore,Hoechst dye can be used for living cell labeling. The fluorescence intensity of Hoechst dye increases with the increase of pH of solution ^[1] .
In Vitro	General Protocol Preparation of Hoechst working solution 1.1 Preparation of the stock solution Dissolve 10 mg of in 5 mL DMSO Note: It is recommended to store the stock solution at 4℃ or -20℃ away from light and avoid repetitive freeze-thaw cycles. 1.2 Preparation of Hoechst working solution Dilute the stock solution in serum-free cell culture medium or PBS to obtain final concentration 10 μg/mL Hoechst working solution. Note: Please adjust the concentration of Hoechst working solution according to the actual situation.

1. Cell staining
- 2.1 Suspension cells 6-well plate
 - a. Centrifuge at 1000 g at 4℃ for 3-5 minutes and then discard the supernatant. Wash twice with PBS, 5 minutes each time. The cell density is 1×10^6 /mL.
 - b. Add 1 mL of working solution, and then incubate at room temperature for 3-10 minutes.
 - c. Centrifuge at 400 g at 4℃ for 3-4 minutes and then discard the supernatant.
 - d. Wash twice with PBS, 5 minutes each time.
 - e. Resuspend cells with serum-free cell culture medium or PBS. Observation by fluorescence microscopy or flow cytometry.
- 2.2 Adherent cells
 - a. Culture adherent cells on sterile coverslips.
 - b. Remove the coverslip from the medium and aspirate excess medium.
 - c. Add 100 μ L of working solution, gently shake it to completely cover the cells, and then incubate at room temperature for 3-10 minutes.
 - d. Wash twice with medium, 5 minutes each time. Observation by fluorescence microscopy or flow cytometry.

Precautions

1. Please adjust the concentration of Hoechst working solution according to the actual situation.
 2. This product is for R&D use only, not for drug, household, or other uses.
 3. For your safety and health, please wear a lab coat and disposable gloves to operate.
- MCE has not independently confirmed the accuracy of these methods. They are for reference only.

REFERENCES

[1]. Chazotte B. Labeling nuclear DNA with hoechst 33342. Cold Spring Harb Protoc. 2011 Jan 1;2011(1):pdb.prot5557.

Caution: Product has not been fully validated for medical applications. For research use only.

Tel: 609-228-6898

Fax: 609-228-5909

E-mail: tech@MedChemExpress.com

Address: 1 Deer Park Dr, Suite Q, Monmouth Junction, NJ 08852, USA