



SZABO SCANDIC

Part of Europa Biosite

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!
See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC Handels GmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7


mail@szabo-scandic.com

www.szabo-scandic.com

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 

Forodesine

Cat. No.:	HY-16210
CAS No.:	209799-67-7
Molecular Formula:	C ₁₁ H ₁₄ N ₄ O ₄
Molecular Weight:	266.25
Target:	Nucleoside Antimetabolite/Analog; Apoptosis
Pathway:	Cell Cycle/DNA Damage; Apoptosis
Storage:	Powder -20°C 3 years 4°C 2 years In solvent -80°C 6 months -20°C 1 month

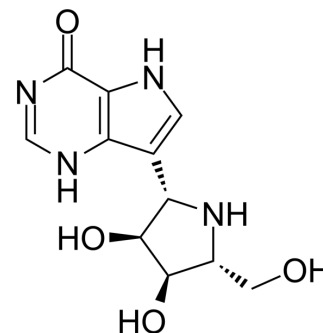

SOLVENT & SOLUBILITY

In Vitro	DMSO : 100 mg/mL (375.59 mM; Need ultrasonic)					
	Preparing Stock Solutions	Solvent Concentration	Mass	1 mg	5 mg	10 mg
		1 mM		3.7559 mL	18.7793 mL	37.5587 mL
		5 mM		0.7512 mL	3.7559 mL	7.5117 mL
		10 mM		0.3756 mL	1.8779 mL	3.7559 mL
Please refer to the solubility information to select the appropriate solvent.						
In Vivo	1. Add each solvent one by one: 10% DMSO >> 40% PEG300 >> 5% Tween-80 >> 45% saline Solubility: ≥ 2.5 mg/mL (9.39 mM); Clear solution					
	2. Add each solvent one by one: 10% DMSO >> 90% (20% SBE-β-CD in saline) Solubility: ≥ 2.5 mg/mL (9.39 mM); Clear solution					
	3. Add each solvent one by one: 10% DMSO >> 90% corn oil Solubility: ≥ 2.5 mg/mL (9.39 mM); Clear solution					

BIOLOGICAL ACTIVITY

Description	Forodesine (BCX-1777) is a highly potent and orally active purine nucleoside phosphorylase (PNP) inhibitor with IC ₅₀ values ranging from 0.48 to 1.57 nM for human, mouse, rat, monkey and dog PNP. Forodesine is a potent human lymphocyte proliferation inhibitor. Forodesine could induce apoptosis in leukemic cells by increasing the dGTP levels ^{[1][2]} .
IC ₅₀ & Target	IC ₅₀ : 1.19 nM (Human PNP), 0.48 nM (Mouse PNP), 1.24 nM (Rat PNP), 0.66 nM (Monkey PNP) and 1.57 nM (Dog PNP) ^[2]
In Vitro	Forodesine (10-30 μM; 24 and 48 hours; RPMI-8226, MOLT-4 and 5T33MM cells) treatment is partially inhibition of

proliferation^[1].

Forodesine (10-30 μ M; 24 and 48 hours; RPMI-8226, MOLT-4 and 5T33MM cells) has no effect on the MM cells at 24 hours, while it could reduce the percentage of living cells in the MOLT-4 cells with 40%^[1].

Forodesine (BCX-1777), in the presence of 2'-deoxyguanosine (dGuo, 3-10 μ M), inhibits human lymphocyte proliferation activated by various agents such as interleukin-2 (IL-2), mixed lymphocyte reaction (MLR) and phytohemagglutinin (PHA) (IC50 values < 0.1-0.38 μ M)^[2].

MCE has not independently confirmed the accuracy of these methods. They are for reference only.

Cell Proliferation Assay^[1]

Cell Line:	Human RPMI-8226, human MOLT-4 (T-ALL) cells, 5T33MM (Multiple myeloma, MM)
Concentration:	10 μ M, 20 μ M, 30 μ M
Incubation Time:	24 and 48 hours
Result:	At the effects at 48 hours, a complete block in proliferation in the MOLT-4 cells and 15% reduction in the 5T33MM cells.

Apoptosis Analysis^[1]

Cell Line:	Human RPMI-8226, human MOLT-4 (T-ALL) cells, 5T33MM (Multiple myeloma, MM)
Concentration:	10 μ M, 20 μ M, 30 μ M
Incubation Time:	24 and 48 hours
Result:	A limited induction of apoptosis.

In Vivo

Forodesine (BCX-1777) has excellent oral bioavailability (63%) in mice^[2].

At a single dose of 10 mg/kg in mice, Forodesine elevates dGuo to approximately 5 μ M^[2].

In the human peripheral blood lymphocyte severe combined immunodeficiency (hu-PBL-SCID) mouse model, Forodesine is effective in prolonging the life span 2-fold or more^[2].

MCE has not independently confirmed the accuracy of these methods. They are for reference only.

CUSTOMER VALIDATION

- Cell. 2019 Oct 31;179(4):864-879.e19.
- Nat Commun. 2017 Aug 14;8(1):241.
- Nucleic Acids Res. 2020 Sep 25;48(17):e101.
- Biosens Bioelectron. 2020 Dec 1;169:112616.
- Int J Biol Sci. 2023; 19(3):772-788.

See more customer validations on www.MedChemExpress.com

REFERENCES

[1]. Bieghs L, et al. The effects of forodesine in murine and human multiple myeloma cells. Adv Hematol. 2010;2010:131895.

[2]. Bantia S, et al. Purine nucleoside phosphorylase inhibitor BCX-1777 (Immuicillin-H)--a novel potent and orally active immunosuppressive agent. Int Immunopharmacol. 2001 Jun;1(6):1199-210.

Caution: Product has not been fully validated for medical applications. For research use only.

Tel: 609-228-6898

Fax: 609-228-5909

E-mail: tech@MedChemExpress.com

Address: 1 Deer Park Dr, Suite Q, Monmouth Junction, NJ 08852, USA