



SZABO SCANDIC

Part of Europa Biosite

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!
See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7


mail@szabo-scandic.com

www.szabo-scandic.com

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 

ML-323

Cat. No.:	HY-17543
CAS No.:	1572414-83-5
Molecular Formula:	C ₂₃ H ₂₄ N ₆
Molecular Weight:	384.48
Target:	Deubiquitinase
Pathway:	Cell Cycle/DNA Damage
Storage:	<div> <div>Powder</div> <div>-20°C 3 years</div> <div>4°C 2 years</div> </div> <div> <div>In solvent</div> <div>-80°C 1 year</div> <div>-20°C 6 months</div> </div>

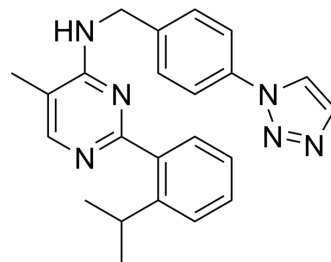

SOLVENT & SOLUBILITY

In Vitro

DMSO : ≥ 49 mg/mL (127.44 mM)
 * "≥" means soluble, but saturation unknown.

	Solvent Concentration	Mass	1 mg	5 mg	10 mg
Preparing Stock Solutions	1 mM		2.6009 mL	13.0046 mL	26.0092 mL
	5 mM		0.5202 mL	2.6009 mL	5.2018 mL
	10 mM		0.2601 mL	1.3005 mL	2.6009 mL

Please refer to the solubility information to select the appropriate solvent.

In Vivo

- Add each solvent one by one: 10% DMSO >> 40% PEG300 >> 5% Tween-80 >> 45% saline
 Solubility: ≥ 2.5 mg/mL (6.50 mM); Clear solution
- Add each solvent one by one: 10% DMSO >> 90% (20% SBE-β-CD in saline)
 Solubility: ≥ 2.5 mg/mL (6.50 mM); Clear solution
- Add each solvent one by one: 10% DMSO >> 90% corn oil
 Solubility: ≥ 2.5 mg/mL (6.50 mM); Clear solution

BIOLOGICAL ACTIVITY

Description

ML-323 is a reversible, potent USP1-UAF1 inhibitor with IC₅₀ of 76 nM in a Ub-Rho assay. The measured inhibition constants of ML-323 for the free enzyme (K_i) is 68 nM.

IC₅₀ & Target

IC₅₀: 76 nM (USP1-UAF1, in Ub-Rho assay)^[1]
 K_i: 68 nM (USP1-UAF1)^[1]

In Vitro

ML-323 (ML323) is a highly potent inhibitor of the USP1-UAF1 deubiquitinase complex with excellent selectivity against human DUBs, deSUMOylase, deneddylase and unrelated proteases. ML-323 is a potent USP1-UAF1 inhibitor with IC₅₀ values of 76 nM in a ubiquitin-rhodamine (Ub-Rho) assay and 174 nM and 820 nM in orthogonal gel-based assays using K63-linked diubiquitin (di-Ub) and monoubiquitinated PCNA (Ub-PCNA) as substrates, respectively. ML-323 probably exerts its inhibitory effect through an allosteric mechanism. The measured inhibition constants of ML-323 for the free enzyme (K_i) and the enzyme-substrate complex (K_i') are 68 nM and 183 nM. Besides, ML-323 potentiates Cisplatin cytotoxicity in non-small cell lung cancer and osteosarcoma cells^[1]. ML-323 (ML323), a probe molecule that displays reversible, nanomolar inhibitory activity and excellent selectivity toward USP1/UAF1. In addition, ML-323 potentiates the cytotoxicity of Cisplatin and increases endogenous monoubiquitination levels of both PCNA and FANCD2, two known cellular targets of USP1/UAF1^[2]. MCE has not independently confirmed the accuracy of these methods. They are for reference only.

PROTOCOL

Kinase Assay ^[1]

For DUB profiling, ML-323 is tested at a single-dose of 10 μ M in duplicate. The DUB activities are monitored using Ub-7-amido-4-methylcoumarin (AMC) as a substrate. The increase in fluorescent signal from free AMC is monitored over time, although only the initial linear portion of slope (signal/min) is used for analysis. The activity of enzyme with no compound is treated as 100%. For protease profiling, ML-323 is tested using threefold serial dilutions starting at 20 μ M against 70 proteases. Proteases are pre-incubated with the compound for 5-15 min before the addition of the appropriate enzyme substrates. The enzyme activities are measured by reading the fluorescent signal from fluorescently labeled peptides^[1]. MCE has not independently confirmed the accuracy of these methods. They are for reference only.

Cell Assay ^[1]

Cell viability is measured by a cell counting kit (CCK) assay using CCK-8 solution. For the colony-forming assay, cells are seeded at a density of 300-500 cells per well in six-well plates and grown overnight. Cells are then treated with ML-323 alone, Cisplatin alone or a combination of Cisplatin and ML-323 (1:1 or 1:4) at the indicated concentrations. Cells treated with an equal volume of DMSO and saline are used as a control. After 48 h of treatment, fresh growth medium is added, and cells are incubated for an additional 5-10 d to allow for colony formation. For UV combination treatment, the cells are treated with ML-323 at the indicated concentrations or an equal volume of DMSO. After 48 h, the medium is removed, and cells are irradiated at 254 nm at the indicated dosage. Fresh growth medium is added, and the cells are incubated for an additional 5-10 d to allow for colony formation. The cells without UV irradiation but treated with ML-323 or an equal volume of DMSO are used as controls and designated as 100%. After the formation of the colonies, cells are fixed with methanol and stained with 0.5% crystal violet. Colonies consisting of >50 cells are scored. The number of colonies is determined from triplicate plates. The dose-response curves are generated using GraphPad Prism and analyzed by using CalcuSyn to calculate the combination index, which is determined for the fraction of cells affected after the addition of fixed ratios of cisplatin and the USP1-UAF1 inhibitor^[1]. MCE has not independently confirmed the accuracy of these methods. They are for reference only.

CUSTOMER VALIDATION

- Nat Commun. 2022 Mar 31;13(1):1700.
- EMBO J. 2015 May 12;34(10):1385-98.
- Int J Biol Sci. 2022 May 1;18(8):3122-3136.
- Cell Death Dis. 2022 Nov 10;13(11):951.
- Cell Chem Biol. 2021 Apr 27;S2451-9456(21)00213-0.

See more customer validations on www.MedChemExpress.com

REFERENCES

[1]. Liang Q, et al. A selective USP1-UAF1 inhibitor links deubiquitination to DNA damage responses. Nat Chem Biol. 2014 Apr;10(4):298-304.

[2]. Dexheimer TS, et al. Discovery of ML323 as a Novel Inhibitor of the USP1/UAF1 Deubiquitinase Complex. National Center for Biotechnology Information; 2010-2012 Oct 23.

Caution: Product has not been fully validated for medical applications. For research use only.

Tel: 609-228-6898

Fax: 609-228-5909

E-mail: tech@MedChemExpress.com

Address: 1 Deer Park Dr, Suite Q, Monmouth Junction, NJ 08852, USA