



# SZABO SCANDIC

Part of Europa Biosite

## Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!  
See the following pages for more information!



### Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

### Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

### SZABO-SCANDIC Handels GmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7


[mail@szabo-scandic.com](mailto:mail@szabo-scandic.com)

[www.szabo-scandic.com](http://www.szabo-scandic.com)

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 

## L-Thyroxine sodium salt pentahydrate

Cat. No.:	HY-18341A
CAS No.:	6106-07-6
Molecular Formula:	C <sub>15</sub> H <sub>20</sub> I <sub>4</sub> NNaO <sub>9</sub>
Molecular Weight:	888.93
Target:	Thyroid Hormone Receptor; Endogenous Metabolite
Pathway:	Vitamin D Related/Nuclear Receptor; Metabolic Enzyme/Protease
Storage:	4°C, sealed storage, away from moisture * In solvent : -80°C, 6 months; -20°C, 1 month (sealed storage, away from moisture)

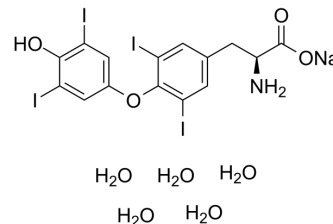

### SOLVENT & SOLUBILITY

In Vitro	DMSO : ≥ 28 mg/mL (31.50 mM) H <sub>2</sub> O : < 0.1 mg/mL (insoluble) * "≥" means soluble, but saturation unknown.					
	Preparing Stock Solutions	<div>Solvent Concentration</div>	Mass	1 mg	5 mg	10 mg
		1 mM		1.1249 mL	5.6247 mL	11.2495 mL
		5 mM		0.2250 mL	1.1249 mL	2.2499 mL
		10 mM		0.1125 mL	0.5625 mL	1.1249 mL
Please refer to the solubility information to select the appropriate solvent.						
In Vivo	1. Add each solvent one by one: 10% DMSO >> 40% PEG300 >> 5% Tween-80 >> 45% saline Solubility: ≥ 2.5 mg/mL (2.81 mM); Clear solution					
	2. Add each solvent one by one: 10% DMSO >> 90% (20% SBE-β-CD in saline) Solubility: 2.5 mg/mL (2.81 mM); Suspended solution; Need ultrasonic					

### BIOLOGICAL ACTIVITY

Description	L-Thyroxine sodium salt pentahydrate (Levothyroxine; T4) is a synthetic hormone for the research of hypothyroidism. DIO enzymes convert biologically active thyroid hormone (Triiodothyronine, T3) from L-Thyroxine (T4) <sup>[1]</sup> .
IC <sub>50</sub> & Target	Human Endogenous Metabolite
In Vivo	Deiodinases (DIOs), which catalyse the conversion of thyroxine (pro-hormone) to the active thyroid hormone, are associated with thyroid stimulating hormone (TSH) levels. DIO1 and DIO2 catalyze activation of thyroid hormone secretion in contrast to DIO3 playing role inactivation of the secretion. Activities of DIO1 and DIO2 play pivotal role in the negative feedback regulation of pituitary TSH secretion <sup>[1]</sup> . L-Thyroxine (T4) and Triiodothyronine (T3) hormones are known to modulate the

expression of ionic channels, pumps and regulatory contractile proteins. Moreover, thyroid hormones have been shown to influence calcium homeostasis and flux responsible for excitation and contractility, with L-Thyroxine and Triiodothyronine modulating its pharmacological control and secretion. In rats fed 12 weeks with the iodine-free diet, a significant decrease in the levels of both Triiodothyronine and L-Thyroxine is observed when compared to the control group fed with standard diet ( $p < 0.001$ ). In the group treated with low doses of L-Thyroxine, an increase in L-Thyroxine levels is observed ( $p = 0.02$ ) while Triiodothyronine levels remain virtually similar to the control group ( $p = 0.19$ ). Rats treated with high doses of L-Thyroxine display a significant increase in both Triiodothyronine and L-Thyroxine circulating concentrations compared to the non-treated hypothyroid group ( $p < 0.001$  and  $p = 0.004$ , respectively) and a significant increase in L-Thyroxine levels when compared to the control values ( $p = 0.03$ )<sup>[2]</sup>.

MCE has not independently confirmed the accuracy of these methods. They are for reference only.

## PROTOCOL

### Animal Administration <sup>[2]</sup>

#### Rats<sup>[2]</sup>

Sprague-Dawley female rats (N=22) are used. Non-pregnant rats are divided into four groups: 1) control, 2) hypothyroidism, 3) hypothyroidism treated with low doses of L-Thyroxine (20 µg/kg/day) and 4) with high doses of L-Thyroxine (100 µg/kg/day). Control rats (group 1) are fed with standard diet, while the intervention rats are fed with iodine-free diet for 12 weeks to induce hypothyroidism (groups 2-4) which is continued for four more weeks to allow screening of hypothyroid status and L-Thyroxine-treatment. Food and water (iodine-free diet) are available ad libitum. The hypothyroid group treated with low (group 3) or high doses of L-Thyroxine (group 4) are injected intraperitoneally every 24 h with respectively 20 µg/kg/day and 100 µg/kg/day. Blood samples are collected for thyroid function screening at week 12 and 16 following the initiation of either the control or iodine-free diet. Hysterectomy is performed under general anesthesia (isoflurane 2%) at the end of the treatment and the two uterine horns are placed in physiological Krebs' solution until isometric tension measurements within no more than 1 h.

MCE has not independently confirmed the accuracy of these methods. They are for reference only.

## CUSTOMER VALIDATION

- Cell Rep Med. 2023 May 24;101061.
- J Drug Deliv Sci Technol. 2023 Sep 28, 105008.
- Stem Cell Rev Rep. 2021 Jun;17(3):999-1013.
- Sci Rep. 2022 Jul 4;12(1):11259.
- Mol Brain. 2021 Jan 27;14(1):25.

See more customer validations on [www.MedChemExpress.com](http://www.MedChemExpress.com)

## REFERENCES

[1]. Arici M, et al. Association between genetic polymorphism and levothyroxine bioavailability in hypothyroid patients. *Endocr J*. 2018 Mar 28;65(3):317-323.

[2]. Corriveau S, et al. Levothyroxine treatment generates an abnormal uterine contractility patterns in an in vitro animal model. *J Clin Transl Endocrinol*. 2015 Sep 9;2(4):144-149.

---

**Caution: Product has not been fully validated for medical applications. For research use only.**

Tel: 609-228-6898

Fax: 609-228-5909

E-mail: [tech@MedChemExpress.com](mailto:tech@MedChemExpress.com)

Address: 1 Deer Park Dr, Suite Q, Monmouth Junction, NJ 08852, USA