



SZABO SCANDIC

Part of Europa Biosite

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!
See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7


mail@szabo-scandic.com

www.szabo-scandic.com

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 

NU2058

Cat. No.:	HY-19316		
CAS No.:	161058-83-9		
Molecular Formula:	C ₁₂ H ₁₇ N ₅ O		
Molecular Weight:	247.3		
Target:	CDK		
Pathway:	Cell Cycle/DNA Damage		
Storage:	Powder	-20°C	3 years
		4°C	2 years
	In solvent	-80°C	2 years
		-20°C	1 year

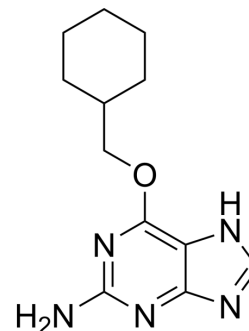

SOLVENT & SOLUBILITY

In Vitro

DMSO : ≥ 32 mg/mL (129.40 mM)

* "≥" means soluble, but saturation unknown.

	Solvent Concentration	Mass	1 mg	5 mg	10 mg
Preparing Stock Solutions	1 mM		4.0437 mL	20.2184 mL	40.4367 mL
	5 mM		0.8087 mL	4.0437 mL	8.0873 mL
	10 mM		0.4044 mL	2.0218 mL	4.0437 mL

Please refer to the solubility information to select the appropriate solvent.

BIOLOGICAL ACTIVITY

Description	NU2058 (O6-(Cyclohexylmethyl)guanine) is a potent, competitive and guanine-based CDK inhibitor with IC ₅₀ s of 17 μM and 26 μM for CDK2 and CDK1. NU2058 has anti-cancer activity.		
IC ₅₀ & Target	CDK1 26 μM (IC ₅₀)	CDK2 17 μM (IC ₅₀)	CDK2 12 μM (Ki)
In Vitro	NU2058 (O6-(Cyclohexylmethyl)guanine) is a competitive inhibitor of CDK2 with respect to ATP (Ki value CDK2, 12 μM) that binds in the ATP binding pocket in a different orientation from other purine-based inhibitors, including olomoucine and roscovitine. NU2058 is the lead compound in a structure-based drug discovery program to develop more potent and selective CDK inhibitors. MCE has not independently confirmed the accuracy of these methods. They are for reference only.		

REFERENCES

Caution: Product has not been fully validated for medical applications. For research use only.

Tel: 609-228-6898

Fax: 609-228-5909

E-mail: tech@MedChemExpress.com

Address: 1 Deer Park Dr, Suite Q, Monmouth Junction, NJ 08852, USA