



SZABO SCANDIC

Part of Europa Biosite

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!
See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC Handels GmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7


mail@szabo-scandic.com

www.szabo-scandic.com

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 

Rf470DL

Cat. No.:	HY-D1689
CAS No.:	2758364-03-1
Molecular Formula:	C ₂₆ H ₃₁ ClN ₄ O ₅ S
Molecular Weight:	547.07
Target:	Bacterial
Pathway:	Anti-infection
Storage:	Please store the product under the recommended conditions in the Certificate of Analysis.

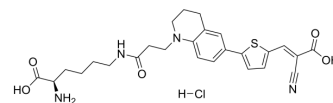

BIOLOGICAL ACTIVITY

Description	Rf470DL is a rotor-fluorogenic D-amino acid (RfDAA). Rf470DL can be used for labeling bacteria (Ex=470 nm, Em=640 nm) ^[1] .
In Vitro	Rf470DL can realize ash-free imaging of bacterial cell walls^[1]. Guidelines (Following is our recommended protocol. This protocol only provides a guideline, and should be modified according to your specific needs)^[1]. 1. RfDAA/FDAA stock solutions are prepared in DMSO at a concentration of 100 mM and stored at -20°C before use. 2.1 For long-pulse labeling of <i>B. subtilis</i> and <i>E. coli</i>, dilute the exponential phase cultures with fresh LB broth containing 1 mM Rf470DL to OD₆₀₀~0.05, and incubate for 1 h. 2.2 Image immediately using a Nikon Ti-E inverted microscopy system without washing and fixation. 3.1 For short-pulse labeling of <i>S. venezuelae</i>, Rf470DL stock solution is added directly to exponential phase cultures to a final concentration of 0.5 mM, followed by incubating at 30 °C with shaking for 15 min. 3.2 Image immediately using a Nikon Ti-E inverted microscopy system without washing and fixation. MCE has not independently confirmed the accuracy of these methods. They are for reference only.

REFERENCES

[1]. Yen-Pang Hsu, et al. Fluorogenic D-amino acids enable real-time monitoring of peptidoglycan biosynthesis and high-throughput transpeptidation assays. *Nat Chem.* 2019 Apr;11(4):335-341.

Caution: Product has not been fully validated for medical applications. For research use only.

Tel: 609-228-6898

Fax: 609-228-5909

E-mail: tech@MedChemExpress.com

Address: 1 Deer Park Dr, Suite Q, Monmouth Junction, NJ 08852, USA