

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten! See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere Liefer- und Versandbedingungen

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

mail@szabo-scandic.com

www.szabo-scandic.com

linkedin.com/company/szaboscandic in





Datasheet for 605-741-125

Goat IgG (H&L) Antibody DyLight™ 488 Conjugated Pre-Adsorbed

Overview

Description:	Donkey Anti-Goat IgG (H&L) Antibody DyLight™ 488 Conjugated (Min X Ch GP Ham Hs Ms Rb 8 Rt Serum Proteins) - 605-741-125			
Item No.:	605-741-125			
Size:	100 μg			
Applications:	Dot Blot, IF, Multiplex			
Reactivity:	Goat			
Host Species:	Donkey			

Product Details

De alcanación de

background:	Anti-Goat igo bylight Antibody generated in donkey detects goat igo. Secreted as part of the
	adaptive immune response by plasma B cells, immunoglobulin G constitutes 75% of serum
	immunoglobulins. Immunoglobulin G binds to viruses, bacteria, as well as fungi and facilitates
	their destruction or neutralization via agglutination (and thereby immobilizing them), activation
	of the compliment cascade, and opsonization for phagocytosis. The whole IgG molecule
	noscesses both the E(c) region, recognized by high-affinity Ec recentor proteins, as well as the E

possesses both the F(c) region, recognized by high-affinity Fc receptor proteins, as well as the F (ab) region possessing the epitope-recognition site. Both heavy and light chains of the antibody

molecule are present.

Synonyms: donkey anti-Goat IgG Antibody DyLight™ 488 Conjugation, donkey anti-Goat IgG DyLight™ 488

Conjugated Antibody

Host Species: Donkey

Specificity: IgG (H&L)

Conjugate: DyLight™ 488

Clonality: Polyclonal

Format: IgG

F/P Ratio: 5.6

Target Details

Reactivity: Goat

www.rockland.com Page 1 of 5



Immunogen:	Goat IgG whole molecule
Purity/Specificity:	This product was prepared from monospecific antiserum by immunoaffinity chromatography using Goat IgG coupled to agarose beads followed by solid phase adsorption(s) to remove any unwanted reactivities. Assay by immunoelectrophoresis resulted in a single precipitin arc against anti-Donkey Serum, Goat IgG and Goat Serum. No reaction was observed against Chicken, Guinea Pig, Hamster, Horse, Mouse, Rabbit and Rat Serum Proteins. This antibody will react with heavy chains of Goat IgG and with light chains of most Goat immunoglobulins.

Application Details

Tested Applications:	Dot Blot			
Suggested Applications:	IF, Multiplex (Based on references)			
Application Note:	Anti-Goat IgG DyLight 488 Antibody has been tested by dot blot. This product is designed for immunofluorescence microscopy, fluorescence based plate assays (FLISA) and fluorescent western blotting. This product is also suitable for multiplex analysis, including multicolor imaging, utilizing various commercial platforms. The emission spectra for this DyLight™ conjugate match the principle output wavelengths of most common fluorescence instrumentation.			
Assay Dilutions: All assays should be optimized by the user. Recommended dilutions (if any) may be listed below.				
FLISA:	>1:20,000			
IF:	>1:5,000			
WB:	>1:10,000			

Formulation

Physical State:	Lyophilized		
Concentration:	1.0 mg/mL by UV absorbance at 280 nm		
Buffer:	0.02 M Potassium Phosphate, 0.15 M Sodium Chloride, pH 7.2		
Preservative:	0.01% (w/v) Sodium Azide		
Stabilizer:	10 mg/mL Bovine Serum Albumin (BSA) - Immunoglobulin and Protease free		
Reconstitution Volume:	100 μL		
Reconstitution Buffer:	Restore with deionized water (or equivalent)		


Shipping & Handling

www.rockland.com Page 2 of 5



Shipping Condition:	Ambient Store vial at 4° C prior to restoration. For extended storage aliquot contents and freeze at -20° C or below. Avoid cycles of freezing and thawing. Centrifuge product if not completely clear after standing at room temperature. This product is stable for several weeks at 4° C as an undiluted liquid. Dilute only prior to immediate use.		
Storage Condition:			
Expiration:	Expiration date is one (1) year from date of receipt.		

Images

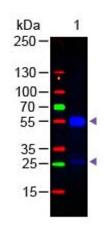

Western Blot


Western Blot of Donkey anti-Goat IgG Antibody DyLight™ 488 Conjugated Pre-Adsorbed.

Lane 1: Goat IgG. Load: 50 ng per lane.

Secondary antibody: Goat IgG (H&L) Antibody DyLight™ 488 Conjugated Pre-Adsorbed at 1:1,000 for 60 min at RT.

Block: MB-070 for 30 min at RT.

Predicted/Observed size: 55 and 28 kDa.

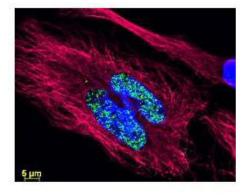

Immunofluorescence Microscopy

DyLight™ dyes can be used for multi-color immunofluorescence microscopy with uniform fluorescence intensity throughout the image. DyLight™ dyes are exceptionally bright and photostable and are optimized for microscopy and microarray detection methods. This image shows anti-Histone detection using a DyLight™ 488 conjugate (green). Anti-Tubulin was detected using a DyLight™ 549 conjugate (red). Nuclei were counter-stained using DAPI (blue). The image was captured using an Axio Imager.Z1 (Zeiss Micro Imaging Inc).

www.rockland.com Page 3 of 5






Diagram

Properties of Dylight™ Conjugates.

Emission	Color	DyLight™ Dye	Ex/Em (nm)	е (M-1 cm-1)	Similar Dyes
Blue		405	400/420	30,000	Alexa™ 405, Cascade Blue
Green		488	493/518	70,000	Alexa™ 488, Cy2®, FITC
Yellow		549	550/568	150,000	Alexa™ 546, Alexa 555, Cy3®,TRITC
Red		649	646/674	250,000	Alexa™ 647, Cy5®
Near Infrared		680	682/715	140,000	Alexa™ 680, Cy5.5®, IRDye™ 700
Infrared		800	770/794	270,000	IRDye™ 800

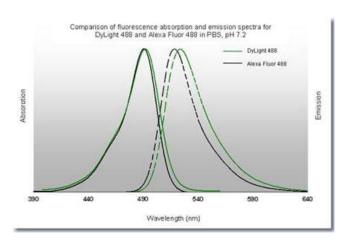

Diagram

References

www.rockland.com Page 4 of 5





- Onoda A et al. A Novel Staining Method for Detection of Brain Perivascular Injuries Induced by Nanoparticle: Periodic Acid-Schiff and Immunohistochemical Double-Staining. Front Toxicol. (2022)
- Iguchi, N et al. Functional constipation induces bladder overactivity associated with upregulations of Htr2 and Trpv2 pathways. *Scientific Reports* (2021)
- Onoda A, Kawasaki T, Tsukiyama K, Takeda K, Umezawa M. Carbon nanoparticles induce endoplasmic reticulum stress around blood vessels with accumulation of misfolded proteins in the developing brain of offspring. *Sci Rep.* (2020)
- Umezawa M et al. Maternal inhalation of carbon black nanoparticles induces neurodevelopmental changes in mouse offspring. Particle and Fibre Toxicology (2018)
- Onoda et al. Dose-dependent induction of astrocyte activation and reactive astrogliosis in mouse brain following maternal exposure to carbon black nanoparticle. Particle and Fibre Toxicology (2017)
- Onoda et al. Perivascular Accumulation of β-Sheet-Rich Proteins in Offspring Brain following Maternal Exposure to Carbon Black Nanoparticles. *Frontiers in Cellular Neuroscience* (2017)
- · Nishino T et al. Distinct distribution of the tensin family in the mouse kidney and small intestine. Exp Anim. (2012)

Disclaimer

This product is for research use only and is not intended for therapeutic or diagnostic applications. Please contact a technical service representative for more information. All products of animal origin manufactured by Rockland Immunochemicals are derived from starting materials of North American origin. Collection was performed in United States Department of Agriculture (USDA) inspected facilities and all materials have been inspected and certified to be free of disease and suitable for exportation. All properties listed are typical characteristics and are not specifications. All suggestions and data are offered in good faith but without guarantee as conditions and methods of use of our products are beyond our control. All claims must be made within 30 days following the date of delivery. The prospective user must determine the suitability of our materials before adopting them on a commercial scale. Suggested uses of our products are not recommendations to use our products in violation of any patent or as a license under any patent of Rockland Immunochemicals, Inc. If you require a commercial license to use this material and do not have one, then return this material, unopened to: Rockland Inc., P.O. BOX 5199, Limerick, Pennsylvania, USA.

www.rockland.com Page 5 of 5