

Produktinformation



Forschungsprodukte & Biochemikalien
Zellkultur & Verbrauchsmaterial
Diagnostik & molekulare Diagnostik
Laborgeräte & Service

Weitere Information auf den folgenden Seiten! See the following pages for more information!



Lieferung & Zahlungsart siehe unsere Liefer- und Versandbedingungen

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien T. +43(0)1 489 3961-0 F. +43(0)1 489 3961-7 <u>mail@szabo-scandic.com</u> www.szabo-scandic.com



www.rockland.com tech@rockland.com +1 484.791.3823

Datasheet for 605-742-002 Goat IgG (H&L) Antibody DyLight™ 549 Conjugated

Overview

Description:	Donkey Anti-Goat IgG (H&L) Antibody DyLight™ 549 Conjugated - 605-742-002		
Item No.:	605-742-002		
Size:	100 µg		
Applications:	Dot Blot, IF, Multiplex		
Reactivity:	Goat		
Host Species:	Donkey		

Product Details

Background:	Anti-Goat IgG DyLight Antibody generated in donkey detects goat IgG. Secreted as part of the adaptive immune response by plasma B cells, immunoglobulin G constitutes 75% of serum immunoglobulins. Immunoglobulin G binds to viruses, bacteria, as well as fungi and facilitates their destruction or neutralization via agglutination (and thereby immobilizing them), activation of the compliment cascade, and opsonization for phagocytosis. The whole IgG molecule possesses both the F(c) region, recognized by high-affinity Fc receptor proteins, as well as the F (ab) region possessing the epitope-recognition site. Both heavy and light chains of the antibody molecule are present.
Synonyms:	donkey anti-Goat IgG Antibody DyLight™ 549 Conjugation, donkey anti-Goat IgG DyLight™ 549 Conjugated Antibody
Host Species:	Donkey
Specificity:	IgG (H&L)
Conjugate:	DyLight™ 549
Clonality:	Polyclonal
Format:	lgG
F/P Ratio:	2.6

Target Details

Reactivity:

Goat



www.rockland.com tech@rockland.com +1 484.791.3823

Immunogen:	Goat IgG, whole molecule		
Purity/Specificity:	This product was prepared from monospecific antiserum by immunoaffinity chromatography using Goat IgG coupled to agarose beads followed by conjugation to fluorochrome and extensive dialysis against the buffer stated above. Assay by immunoelectrophoresis resulted in a single precipitin arc against anti-Donkey Serum, Goat IgG and Goat Serum. This antibody will react with heavy chains of Goat IgG and with light chains of most Goat immunoglobulins.		

Application Details

Tested Applications:	Dot Blot			
Suggested Applications:	IF, Multiplex (Based on references) Anti-Goat IgG DyLight 549 Antibody has been tested by dot blot. This product is designed for immunofluorescence microscopy, fluorescence based plate assays (FLISA) and fluorescent western blotting. This product is also suitable for multiplex analysis, including multicolor imaging, utilizing various commercial platforms. The emission spectra for this DyLight [™] conjugate match the principle output wavelengths of most common fluorescence instrumentation.			
Application Note:				
Assay Dilutions:	All assays should be optimized by the user. Recommended dilutions (if any) may be listed below.			
FLISA:	>1:20,000			
IF:	>1:5,000			
WB:	>1:10,000			

Formulation

Physical State:	Lyophilized		
Concentration:	1.0 mg/mL by UV absorbance at 280 nm		
Buffer:	0.02 M Potassium Phosphate, 0.15 M Sodium Chloride, pH 7.2		
Preservative:	0.01% (w/v) Sodium Azide		
Stabilizer:	10 mg/mL Bovine Serum Albumin (BSA) - Immunoglobulin and Protease free		
Reconstitution Volume:	100 μL		
Reconstitution Buffer:	Restore with deionized water (or equivalent)		

Shipping & Handling

Shipping Condition: Ambient



www.rockland.com tech@rockland.com +1 484.791.3823

Storage Condition:	Store vial at 4° C prior to restoration. For extended storage aliquot contents and freeze at -20° C or below. Avoid cycles of freezing and thawing. Centrifuge product if not completely clear after standing at room temperature. This product is stable for several weeks at 4° C as an undiluted liquid. Dilute only prior to immediate use.
Expiration:	Expiration date is one (1) year from date of receipt.

Images

Diagram

Properties of DyLight[™] Conjugates.

Emission	Color	DyLight™ Dye	Ex/Em (nm)	е (M ⁻¹ cm ⁻¹)	Similar Dyes
Blue		405	400/420	30,000	Alexa™ 405, Cascade Blue
Green		488	493/518	70,000	Alexa™ 488, Cy2®, FITC
Yellow		549	550/568	150,000	Alexa™ 546, Alexa 555, Cy3®,TRITC
Red		649	646/674	250,000	Alexa™ 647, Cy5®
Near Infrared		680	682/715	140,000	<u>Alexa™ 680, Cy5.5®, IRDye™ 700</u>
Infrared		800	770/794	270,000	IRDye™ 800

Tupniu 25 ng 1 ng 6 ng 1 ng 1 ng 25 ng 12 ng 5 ng 12 ng 5 ng 1 ng 1 ng 25 ng 12 ng 5 ng 1 ng 1 ng 25 ng 1 ng 1 ng 25 ng 1 ng 1 ng 1 2 ng 1 2


Western Blot

DyLight[™] dyes can be used for two-color western blot detection with low background and high signal. Anti-tubulin was detected using a DyLight[™] 549 conjugate. Anti-TNFa was detected using a DyLight[™] 649 conjugate. The image was captured using the Typhoon[™] 9410 Imaging System.



Order online now!

www.rockland.com tech@rockland.com +1 484.791.3823

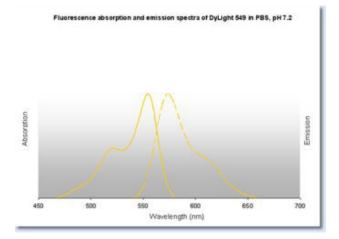

Diagram

DyLight[™] 549 Fluorescence Spectra.

References

- Onoda et al. Perivascular Accumulation of β-Sheet-Rich Proteins in Offspring Brain following Maternal Exposure to Carbon Black Nanoparticles. *Frontiers in Cellular Neuroscience* (2017)
- Nishino T et al. Distinct distribution of the tensin family in the mouse kidney and small intestine. *Exp Anim.* (2012)

Disclaimer

This product is for research use only and is not intended for therapeutic or diagnostic applications. Please contact a technical service representative for more information. All products of animal origin manufactured by Rockland Immunochemicals are derived from starting materials of North American origin. Collection was performed in United States Department of Agriculture (USDA) inspected facilities and all materials have been inspected and certified to be free of disease and suitable for exportation. All properties listed are typical characteristics and are not specifications. All suggestions and data are offered in good faith but without guarantee as conditions and methods of use of our products are beyond our control. All claims must be made within 30 days following the date of delivery. The prospective user must determine the suitability of our materials before adopting them on a commercial scale. Suggested uses of our products are not recommendations to use our products in violation of any patent or as a license under any patent of Rockland Immunochemicals, Inc. If you require a commercial license to use this material and do not have one, then return this material, unopened to: Rockland Inc., P.O. BOX 5199, Limerick, Pennsylvania, USA.