

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten! See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere Liefer- und Versandbedingungen

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

mail@szabo-scandic.com

www.szabo-scandic.com

linkedin.com/company/szaboscandic in





Datasheet for 609-442-002

Human IgG (H&L) Antibody Dylight™ 549 Conjugated

Overview

Description:	Rabbit Anti-Human IgG (H&L) Antibody DyLight™ 549 Conjugated - 609-442-002		
Item No.:	609-442-002		
Size:	100 μg		
Applications:	Dot Blot, WB		
Reactivity:	Human		
Host Species:	Rabbit		

Product Details

Background: Anti-Human IgG (H&L) DyLight 549 generated in rabbit detects hu	uman Immunoglobulin G (IgG),
--	------------------------------

both heavy and light chains of the antibody molecule are present. It is a protein complex composed of four peptide chains — two identical heavy chains and two identical light chains arranged in a Y-shape typical of antibody monomers. Each IgG has two antigen binding sites. Representing approximately 75% of serum immunoglobulins in humans, IgG is the most abundant antibody isotype found in the circulation. IgG molecules are synthesized and secreted by plasma B cells. Secondary Antibodies are available in a variety of formats and conjugate types. When choosing a secondary antibody product, consideration must be given to species and immunoglobulin specificity, conjugate type, fragment and chain specificity, level of cross-reactivity, and host-species source and fragment composition.

Synonyms:	rabbit anti-Human IgG Antibody DyLight™549 Conjugation, rabbit anti-Human IgG DyLight™ 549
	Conjugated Antibody

	Conjugated Antibody
Host Species:	Rabbit
Specificity:	IgG (H&L)
Conjugate:	DyLight™ 549
Clonality:	Polyclonal
Format:	IgG
F/P Ratio:	1.2

Target Details

www.rockland.com Page 1 of 4



Reactivity:	Human		
Immunogen:	Human IgG, whole molecule		
Purity/Specificity:	This product was prepared from monospecific antiserum by immunoaffinity chromatography using Human IgG coupled to agarose beads followed by conjugation to fluorochrome and extensive dialysis against the buffer stated above. Assay by immunoelectrophoresis resulted in a single precipitin arc against anti-Rabbit Serum, Human IgG and Human Serum. This antibody will react with heavy chains of Human IgG and with light chains of most Human immunoglobulins.		

Application Details

Tested Applications:	Dot Blot, WB
Application Note:	Anti-Human IgG (H&L) DyLight 549 has been tested by dot blot and western blot and is designed for immunofluorescence microscopy, fluorescence based plate assays (FLISA) and fluorescent western blotting. This product is also suitable for multiplex analysis, including multicolor imaging, utilizing various commercial platforms. The emission spectra for this DyLight™ conjugate match the principle output wavelengths of most common fluorescence instrumentation.
Assay Dilutions:	All assays should be optimized by the user. Recommended dilutions (if any) may be listed below.
FLISA:	>1:20,000
IF:	>1:5,000
WB:	>1:10,000

Formulation

Concentration: 1.0 mg/mL by UV absorbance at 280 nm Buffer: 0.02 M Potassium Phosphate, 0.15 M Sodium Chloride, pH 7.2 Preservative: 0.01% (w/v) Sodium Azide Stabilizer: 10 mg/mL Bovine Serum Albumin (BSA) - Immunoglobulin and Protease free Reconstitution Volume: 100 μL Reconstitution Buffer: Restore with deionized water (or equivalent)	Physical State:	Lyophilized	
Preservative: 0.01% (w/v) Sodium Azide Stabilizer: 10 mg/mL Bovine Serum Albumin (BSA) - Immunoglobulin and Protease free Reconstitution Volume: 100 μL	Concentration:	1.0 mg/mL by UV absorbance at 280 nm	
Stabilizer: 10 mg/mL Bovine Serum Albumin (BSA) - Immunoglobulin and Protease free Reconstitution Volume: 100 μL	Buffer:	0.02 M Potassium Phosphate, 0.15 M Sodium Chloride, pH 7.2	
Reconstitution Volume: 100 μL	Preservative:	0.01% (w/v) Sodium Azide	
	Stabilizer:	10 mg/mL Bovine Serum Albumin (BSA) - Immunoglobulin and Protease free	
Reconstitution Buffer: Restore with deionized water (or equivalent)	Reconstitution Volume:	100 μL	
	Reconstitution Buffer:	Restore with deionized water (or equivalent)	


Shipping & Handling

www.rockland.com Page 2 of 4



Shipping Condition:	Ambient		
Storage Condition:	Store vial at 4° C prior to restoration. For extended storage aliquot contents and freeze at -20° C or below. Avoid cycles of freezing and thawing. Centrifuge product if not completely clear after standing at room temperature. This product is stable for several weeks at 4° C as an undiluted liquid. Dilute only prior to immediate use.		
Expiration:	Expiration date is one (1) year from date of receipt.		

Images

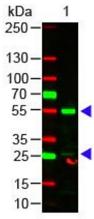

		NDa
		250 —
		130-
		100-
		70 —
4	-	55 —
		35 —
4	•	25-
		15-
		10-

Emission	Color	DyLight™ Dye	Ex/Em (nm)	е (M ⁻¹ cm ⁻¹)	Similar Dyes
Blue		405	400/420	30,000	Alexa™ 405, Cascade Blue
Green	H W	488	493/518	70,000	Alexa™ 488, Cy2®, FITC
Yellow		549	550/568	150,000	Alexa™ 546, Alexa 555, Cy3®,TRITC
Red		649	646/674	250,000	Alexa™ 647, Cy5®
Near Infrared		680	682/715	140,000	Alexa™ 680, Cy5.5®, IRDye™ 700
Infrared		800	770/794	270,000	IRDye™ 800

Western Blot

Western Blot of Rabbit anti-Human IgG (H&L) Antibody DyLight™ 549 Conjugated.

Lane 1: Human IgG. Load: 50 ng per lane. Primary Antibody: None.

Secondary antibody: DyLight™ 549 rabbit secondary

antibody at 1:1,000 for 60 min at RT. Block: MB-070 for 30 min at RT.

Predicted/Observed size: 55 and 28 kDa for Rabbit IgG.


Other Band(s): None.

Diagram

Properties of DyLight™ Fluorescent Dyes.


www.rockland.com Page 3 of 4



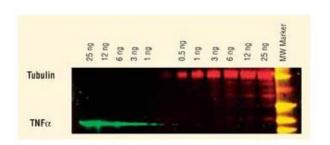

Western Blot

DyLight[™] dyes can be used for two-color Western Blot detection with low background and high signal. Anti-tubulin was detected using a DyLight[™] 549 conjugate. Anti-TNFa was detected using a DyLight[™] 649 conjugate. The image was captured using the Typhoon[™] 9410 Imaging System.

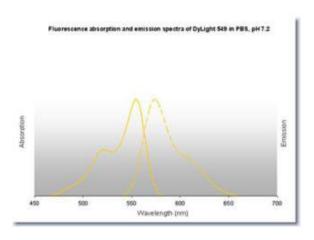

Diagram

DyLight™ 549 Fluorescence Spectra.

Disclaimer

This product is for research use only and is not intended for therapeutic or diagnostic applications. Please contact a technical service representative for more information. All products of animal origin manufactured by Rockland Immunochemicals are derived from starting materials of North American origin. Collection was performed in United States Department of Agriculture (USDA) inspected facilities and all materials have been inspected and certified to be free of disease and suitable for exportation. All properties listed are typical characteristics and are not specifications. All suggestions and data are offered in good faith but without guarantee as conditions and methods of use of our products are beyond our control. All claims must be made within 30 days following the date of delivery. The prospective user must determine the suitability of our materials before adopting them on a commercial scale. Suggested uses of our products are not recommendations to use our products in violation of any patent or as a license under any patent of Rockland Immunochemicals, Inc. If you require a commercial license to use this material and do not have one, then return this material, unopened to: Rockland Inc., P.O. BOX 5199, Limerick, Pennsylvania, USA.

www.rockland.com Page 4 of 4