

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten! See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere Liefer- und Versandbedingungen

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

mail@szabo-scandic.com

www.szabo-scandic.com



www.rockland.com tech@rockland.com +1 484.791.3823

Datasheet for 610-445-042

Mouse IgG2b (Gamma 2b chain) Antibody DyLight™ 800 Conjugated

Overview

Description:	Rabbit Anti-Mouse IgG2b (Gamma 2b chain) Antibody DyLight™ 800 Conjugated - 610-44			
Item No.:	610-445-042			
Size:	100 μg			
Reactivity:	Mouse			
Host Species:	Rabbit			

Product Details

Background: Anti-Mouse IgG2b DyLight 800 Antibody generated in rabbit detects reactivity to Mouse IgG2b

(Gamma 2b chain). Secreted as part of the adaptive immune response by plasma B cells, immunoglobulin G constitutes 75% of serum immunoglobulins. IgG2, the second largest of IgG isotypes, comprises almost 25% of IgG and has a low affinity for binding to the Fc receptor of phagocytic cells. Secondary Antibodies are available in a variety of formats and conjugate types. When choosing a secondary antibody product, consideration must be given to species and immunoglobulin specificity, conjugate type, fragment and chain specificity, level of cross-

reactivity, and host-species source and fragment composition.

Synonyms: Rabbit Anti-Mouse IgG2b (Gamma 2b chain) Antibody DyLight 800™ Conjugated, Rabbit Anti

Mouse IgG2b DyLight 800™ Conjugated Antibody

Host Species: Rabbit

Specificity: IgG2b

Conjugate: DyLight™ 800

Clonality: Polyclonal

Format: IgG

F/P Ratio: 2.0

Target Details

Reactivity: Mouse

Immunogen: Mouse IgG2b heavy chain

www.rockland.com Page 1 of 4





>1:5,000

>1:10,000

www.rockland.com tech@rockland.com +1 484.791.3823

Purity/Specificity:

This product was prepared from monospecific antiserum by immunoaffinity chromatography using antigens coupled to agarose beads followed by solid phase adsorption(s) to remove any unwanted reactivities. Assay by immunoelectrophoresis resulted in a single precipitin arc against anti-Rabbit Serum, Mouse IgG and Mouse Serum. Minimal cross reactivity was noted against other Mouse immunoglobulin classes or subclasses.

Application Details

Application Note:	This product is designed for immunofluorescence microscopy, fluorescence based plate assays (FLISA) and fluorescent western blotting. This product is also suitable for multiplex analysis, including multicolor imaging, utilizing various commercial platforms. The emission spectra for this DyLight™ conjugate match the principle output wavelengths of most common fluorescence instrumentation.		
Assay Dilutions:	All assays should be optimized by the user. Recommended dilutions (if any) may be listed below.		
FLISA:	>1:20,000		

Formulation

IF:

WB:

Physical State:	Lyophilized			
Concentration:	1.0 mg/mL by UV absorbance at 280 nm			
Buffer:	0.02 M Potassium Phosphate, 0.15 M Sodium Chloride, pH 7.2			
Preservative:	0.01% (w/v) Sodium Azide			
Stabilizer:	10 mg/mL Bovine Serum Albumin (BSA) - Immunoglobulin and Protease free			
Reconstitution Volume:	100 μL			
Reconstitution Buffer:	Restore with deionized water (or equivalent)			

Shipping & Handling

Shipping Condition:	Ambient
Storage Condition:	Store vial at 4° C prior to restoration. For extended storage aliquot contents and freeze at -20° C or below. Avoid cycles of freezing and thawing. Centrifuge product if not completely clear after standing at room temperature. This product is stable for several weeks at 4° C as an undiluted liquid. Dilute only prior to immediate use.


www.rockland.com Page 2 of 4



www.rockland.com tech@rockland.com +1 484.791.3823

Expiration: Expiration date is one (1) year from date of receipt.

Images

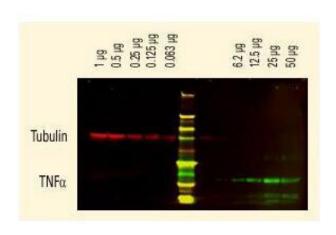

Western Blot

DyLight™ dyes can be used for two-color Western Blot detection with low background and high signal. Anti-tubulin was detected using a DyLight™ 680 conjugate. Anti-TNFa was detected using a DyLight™ 800 conjugate. The image was captured using the Odyssey® Infrared Imaging System developed by LI-COR.

Diagram

Properties of DyLight™ Conjugates.

Emission	Color	DyLight™ Dye	Ex/Em (nm)	е (M ⁻¹ cm ⁻¹)	Similar Dyes
Blue		405	400/420	30,000	Alexa™ 405, Cascade Blue
Green		488	493/518	70,000	Alexa™ 488, Cy2®, FITC
Yellow		549	550/568	150,000	Alexa™ 546, Alexa 555, Cy3®,TRITC
Red		649	646/674	250,000	Alexa™ 647, Cy5®
Near Infrared		680	682/715	140,000	Alexa™ 680, Cy5.5®, IRDye™ 700
Infrared		800	770/794	270,000	IRDye™ 800

www.rockland.com Page 3 of 4



www.rockland.com tech@rockland.com +1 484.791.3823

Fluerescence absorption and emission spectra of DyLight 800 in PBS, pH 7.2

Wavelength (nm)

Diagram

Disclaimer

This product is for research use only and is not intended for therapeutic or diagnostic applications. Please contact a technical service representative for more information. All products of animal origin manufactured by Rockland Immunochemicals are derived from starting materials of North American origin. Collection was performed in United States Department of Agriculture (USDA) inspected facilities and all materials have been inspected and certified to be free of disease and suitable for exportation. All properties listed are typical characteristics and are not specifications. All suggestions and data are offered in good faith but without guarantee as conditions and methods of use of our products are beyond our control. All claims must be made within 30 days following the date of delivery. The prospective user must determine the suitability of our materials before adopting them on a commercial scale. Suggested uses of our products are not recommendations to use our products in violation of any patent or as a license under any patent of Rockland Immunochemicals, Inc. If you require a commercial license to use this material and do not have one, then return this material, unopened to: Rockland Inc., P.O. BOX 5199, Limerick, Pennsylvania, USA.

www.rockland.com Page 4 of 4