

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten! See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere Liefer- und Versandbedingungen

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

mail@szabo-scandic.com

www.szabo-scandic.com



Datasheet for 612-145-002

Rat IgG (H&L) Antibody Dylight™ 800 Conjugated

Overview

Description:	Goat Anti-Rat IgG (H&L) Antibody DyLight™ 800 Conjugated - 612-145-002			
Item No.:	612-145-002			
Size:	100 μg			
Applications:	WB			
Reactivity:	Rat			
Host Species:	Goat			

Product Details

Background:	This product is designed for immunofluorescence microscopy, fluorescence based plate assays (FLISA) and fluorescent western blotting. This product is also suitable for multiplex analysis, including multicolor imaging, utilizing various commercial platforms.			
Synonyms:	Goat Anti-Rat IgG DyLight 800™ Conjugated Antibody, Goat Anti-Rat IgG Antibody DyLight 800™ Conjugation			
Host Species:	Goat			
Specificity:	IgG (H&L)			
Conjugate:	DyLight™ 800			
Clonality:	Polyclonal			
Format:	IgG			
F/P Ratio:	2.0			

Target Details

Reactivity:	Rat
Immunogen:	Rat IgG, whole molecule

www.rockland.com Page 1 of 5





Purity/Specificity: This product was

This product was prepared from monospecific antiserum by immunoaffinity chromatography using Rat IgG coupled to agarose beads followed by conjugation to fluorochrome and extensive dialysis against the buffer stated above. Assay by immunoelectrophoresis resulted in a single precipitin arc against anti-Goat Serum, Rat IgG and Rat Serum. This antibody will react with heavy chains of Rat IgG and with light chains of most Rat immunoglobulins.

Application Details

Tested Applications:	WB		
Application Note:	Anti-Rat IgG DyLight™ 800 has been tested by western blot and is designed for immunofluorescence microscopy, fluorescence based plate assays (FLISA) and fluorescent western blotting. The emission spectra for this DyLight™ conjugate match the principle output wavelengths of most common fluorescence instrumentation.		
Assay Dilutions:	All assays should be optimized by the user. Recommended dilutions (if any) may be listed below.		
FLISA:	>1:20,000		
IF:	>1:5,000		
WB:	>1:10,000		

Formulation

Physical State:	Lyophilized		
Concentration:	1.0 mg/mL by UV absorbance at 280 nm		
Buffer:	0.02 M Potassium Phosphate, 0.15 M Sodium Chloride, pH 7.2		
Preservative:	0.01% (w/v) Sodium Azide		
Stabilizer:	10 mg/mL Bovine Serum Albumin (BSA) - Immunoglobulin and Protease free		
Reconstitution Volume:	100 μL		
Reconstitution Buffer:	Restore with deionized water (or equivalent)		

Shipping & Handling


Shipping Condition:	Ambient
Storage Condition:	Store vial at 4° C prior to restoration. For extended storage aliquot contents and freeze at -20° C or below. Avoid cycles of freezing and thawing. Centrifuge product if not completely clear after standing at room temperature. This product is stable for several weeks at 4° C as an undiluted liquid. Dilute only prior to immediate use.

www.rockland.com Page 2 of 5




Expiration: Expiration date is one (1) year from date of receipt.

Images

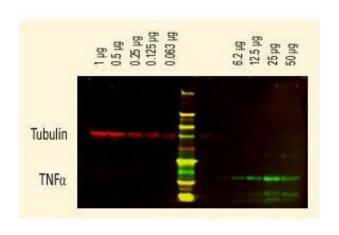

Western Blot

DyLight™ dyes can be used for two-color western blot detection with low background and high signal. Anti-tubulin was detected using a DyLight™ 680 conjugate. Anti-TNFa was detected using a DyLight™ 800 conjugate. The image was captured using the Odyssey® Infrared Imaging System developed by LI-COR.

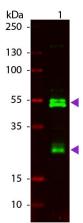


Western Blot

Western Blot of Goat anti-Rat IgG DyLight™ 800 Conjugated Secondary Antibody. Lane 1: Rat IgG. Lane 2: None. Load: 50 ng per lane. Primary antibody: None. Secondary antibody: DyLight™ 800 goat secondary antibody at 1:5,000 for 60 min at RT. Block: MB-070 for 30 min at RT. Predicted/Observed size: 25 & 55 kDa, 25 & 55 kDa for Rat IgG. Other band(s): None.

www.rockland.com Page 3 of 5



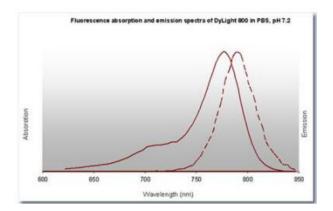

Diagram

DyLight™ 800 Fluorescence Spectra.

Diagram

Properties of DyLight™ Fluorescent Dyes.

Emission	Color	DyLight™ Dye	Ex/Em (nm)	е (M ⁻¹ cm ⁻¹)	Similar Dyes
Blue		405	400/420	30,000	Alexa™ 405, Cascade Blue
Green	7	488	493/518	70,000	Alexa™ 488, Cy2®, FITC
Yellow		549	550/568	150,000	Alexa™ 546, Alexa 555, Cy3®,TRITC
Red		649	646/674	250,000	Alexa™ 647, Cy5®
Near Infrared		680	682/715	140,000	<u>Alexa™ 680, Cy5.5®, IRDye™ 700</u>
Infrared		800	770/794	270,000	IRDye™ 800

References

- Szewczyk B et al. FUS ALS neurons activate major stress pathways and reduce translation as an early protective mechanism against neurodegeneration. *Cell Rep.* (2023)
- Voß Y et al. Malaria parasite centrins can assemble by Ca2+-inducible condensation. PLoS P et al. (2023)
- Balic JJ et al. STAT3 serine phosphorylation is required for TLR4 metabolic reprogramming and IL-1β expression. *Nat Commun.* (2020)
- Bacman SR et al. MitoTALEN reduces mutant mtDNA load and restores tRNA Ala levels in a mouse model of heteroplasmic mtDNA mutation. *Nat Med.* (2018)
- Bacman, SR et al. MitoTALEN reduces mutant mtDNA load and restores tRNAAla levels in a mouse model of heteroplasmic mtDNA mutation. *Nature Medicine* (2018)

www.rockland.com Page 4 of 5





Disclaimer

This product is for research use only and is not intended for therapeutic or diagnostic applications. Please contact a technical service representative for more information. All products of animal origin manufactured by Rockland Immunochemicals are derived from starting materials of North American origin. Collection was performed in United States Department of Agriculture (USDA) inspected facilities and all materials have been inspected and certified to be free of disease and suitable for exportation. All properties listed are typical characteristics and are not specifications. All suggestions and data are offered in good faith but without guarantee as conditions and methods of use of our products are beyond our control. All claims must be made within 30 days following the date of delivery. The prospective user must determine the suitability of our materials before adopting them on a commercial scale. Suggested uses of our products are not recommendations to use our products in violation of any patent or as a license under any patent of Rockland Immunochemicals, Inc. If you require a commercial license to use this material and do not have one, then return this material, unopened to: Rockland Inc., P.O. BOX 5199, Limerick, Pennsylvania, USA.

www.rockland.com Page 5 of 5