

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten! See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere Liefer- und Versandbedingungen

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

mail@szabo-scandic.com

www.szabo-scandic.com

linkedin.com/company/szaboscandic in





Datasheet for WM3772F-01-0010 WM3772F Viable Cells

Overview

| Description: | WM3772F Viable Cells - WM3772F-01-0010 |
|--------------|--|
| Item No.: | WM3772F-01-0010 |
| Size: | 10 x 1 million cells |
| Origin: | Human |

Product Details

| Background: | WM3772F is a primary human metastatic uveal melanoma cell line displaying clumpy pigmented |
|-------------|---|
| | morphology. These cells require 5% FBS to adhere to culture flask. This cell line was established |

morphology. These cells require 5% FBS to adhere to culture flask. This cell line was established from a metastases arising from uveal melanoma of a patient. Uveal melanoma is a rare type of cancer of an eye and can arise in the anterior (iris) or the posterior (ciliary body or choroid) uveal tract. Most uveal tract melanomas originate in the choroid. The ciliary body is less commonly a site of origin, and the iris is the least common. This cell line expresses wild type BRAF, N-RAS, c-KIT, and CDK4 genes. WM3772F cells match with WM3618F. These cells produce

xenograft tumors when injected into immunocompromised mice.

| Synonyms: | Melanoma, patient derived tumor, tumor models, skin cancer, xenograft |
|--------------------|---|
| Species of Origin: | Human |

Target Details

Purity/Specificity: Cells are sterile, validated by short tandem repeat profiling, and are tested as negative for

mycoplasma. It is recommended that cell lines are tested for mycoplasma contamination and short tandem repeat (STR) profiling every 10 passages or each time a frozen seed stock is made.

See cell culture protocol for additional details.

Relevant Links: • Cell Line EULA

Melanoma Cell Culture Protocol

Application Details

Application Note: The key applications of these cell lines include genetic studies, xenograft production, drug

testing, and drug target discovery. These cell line models can be used in various biological

www.rockland.com Page 1 of 4





| | | | 11 | -1 | |
|-----------------------|----------------|-----------------|----------|-----------|-----------|
| assavs, and for ident | itving critica | i target genes. | and cell | signaling | nathways. |

| Assay Dilutions: | All assays should be optimized by the user. Recommended dilutions (if any) may be | | |
|-------------------------|---|--|--|
| | listed below. | | |

Cell Line Data

| Cell Line: | Human Melanoma |
|-----------------|---|
| Product Type: | Viable Cells |
| Morphology: | clumpy pigmented |
| Cell Viability: | Yes |
| Stage: | PRIMARY, Metastasis OF UVEAL |
| BRAF: | WT |
| CDK4: | WT |
| C-Kit: | WT |
| N-RAS: | WT |
| PTEN: | ND |
| Paired: | Yes |
| Medium: | Tumor Specialized Media with 5% HI-FBS |
| Sub-culture: | Cells should be maintained between $30-95\%$ confluence in tumor specialized medium with 5% FBS; split cultures 1:4 every 1 week using 0.25% trypsin/EDTA. |
| Incubation: | 36°C with 5% CO2 |

Formulation

| Concentration: 1.0 million cells/mL Count By Hemocytometer | |
|--|--|
| Concentration: 1.0 million cells/mL Count By Hemocytometer | |
| Buffer: None | |
| Preservative: None | |
| Stabilizer: None | |

Shipping & Handling

Shipping Condition: Dry Ice

www.rockland.com Page 2 of 4






Storage Condition: Cells are frozen with 90% FBS/10% DMSO solution at about 10x10^6 cells/ml. Store vial in liquid

nitrogen upon arrival.

Expiration: Expiration date is two (2) years from date of receipt.

Images

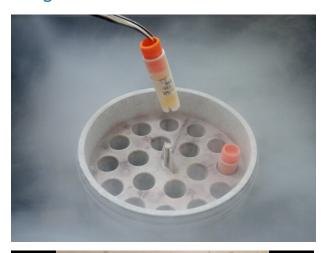

Flask

Human melanoma tumor cells with known gene mutations, disease stage, STR, and RPPA profiling



Viable cell growth

Established WM3772F viable cell growth in culture using appropriate Tumor Specialized Media with 5%FBS.

Disclaimer

www.rockland.com Page 3 of 4





No test method can provide total assurance that the hepatitis B virus, hepatitis C virus, human immunodeficiency virus, or any other infectious agents are absent. Thus, all blood products, including purified proteins derived from human blood sources, should be handled at Biosafety Level 2 as recommended by the CDC\NIH manual entitled Biosafety in Microbiological and Biomedical Laboratories for potentially infectious human serum, blood specimens or proteins derived from same. Source material for the human blood product supplied to your facility has been tested for the detection of HIV antibody, Hepatitis B surface antigen, antibody to Hepatitis C, HIV 1 antigen(s), antibody to HTLV - I/II, and syphilis by FDA guidelines. All units were found to be non-reactive/negative for these tests. All human blood source material is collected in FDA licensed centers and is tested with FDA approved test kits.

Cell Line Limited Use License Required. THIS PRODUCT IS SUBJECT TO AN END-USER LICENSE AGREEMENT (EULA). BY ACCEPTING THIS PRODUCT, RECIPIENT AGREES TO BE BOUND BY THE TERMS OF USE SET FORTH BELOW and SET FORTH IN THE EULA. THIS PRODUCT IS FOR IN VITRO RESEARCH USE ONLY. THERAPEUTIC, DIAGNOSTIC, OR VETERINARY USE IS PROHIBITED. This product may not be resold or transferred by the recipient and may be used only by the recipient, in the recipient's facility and only for research use and other uses specifically permitted by the EULA. No other commercial use is allowed. "Commercial Use" means any and all uses of this product by recipient or others for monetary or other consideration, including providing services, supplying information or data to unaffiliated third parties, and resale or transfer of this product for any use. Recipient has no right to modify, derivatize, genetically engineer or otherwise create variations of this product or associated cells or cell lines. ROCKLAND AND WISTAR MAKE NO REPRESENTATIONS AND EXTEND NO WARRANTIES OF ANY KIND, EITHER EXPRESSED OR IMPLIED, INCLUDING AS TO MERCHANTABILITY OR FITNESS FOR A PARTICULAR PURPOSE, OR THAT THE USE OF THE PRODUCTS WILL NOT INFRINGE ANY PATENT, COPYRIGHT, TRADEMARK, OR OTHER PROPRIETARY RIGHTS. The terms set forth herein and in the EULA shall be governed by the laws of the Commonwealth of Pennsylvania, USA. To obtain a COMMERCIAL USE license for this product, please contact Rockland Immunochemicals, Inc. Please contact a technical service representative for more information. All properties listed are typical characteristics and are not specifications. All suggestions and data are offered in good faith but without guarantee as conditions and methods of use of our products are beyond our control. All claims must be made within 30 days following the date of delivery. The prospective user must determine the suitability of our materials before adopting them on a commercial scale. Suggested uses of our products are not recommendations to use our products in violation of any patent or as a license under any patent of Rockland Immunochemicals, Inc. If you require a commercial license to use this material and do not have one, then return this material, unopened to: Rockland Inc., P.O. BOX 5199, Limerick, Pennsylvania, USA.

www.rockland.com Page 4 of 4