



**SZABO  
SCANDIC**

Part of Europa Biosite

## Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!  
See the following pages for more information!



### Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

### Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

### SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

[mail@szabo-scandic.com](mailto:mail@szabo-scandic.com)

[www.szabo-scandic.com](http://www.szabo-scandic.com)

[linkedin.com/company/szaboscandic](http://linkedin.com/company/szaboscandic)



# PKC $\gamma$ (h): 293T Lysate: sc-116200

## BACKGROUND

Members of the protein kinase C (PKC) family play a key regulatory role in a variety of cellular functions including cell growth and differentiation, gene expression, hormone secretion and membrane function. PKCs were originally identified as serine/threonine protein kinases whose activity was dependent on calcium and phospholipids. Diacylglycerols (DAG) and tumor promoting phorbol esters bind to and activate PKC. PKCs can be subdivided into many different isoforms ( $\alpha$ ,  $\beta I$ ,  $\beta II$ ,  $\gamma$ ,  $\delta$ ,  $\epsilon$ ,  $\zeta$ ,  $\eta$ ,  $\theta$ ,  $\iota$ ,  $\lambda$ ,  $\mu$ , and  $\nu$ ). Patterns of expression for each PKC isoform differ among tissues and PKC family members exhibit clear differences in their cofactor dependencies. For instance, the kinase activities of PKC  $\delta$  and  $\epsilon$  are independent of  $Ca^{2+}$ . On the other hand, most of the other PKC members possess phorbol ester-binding activities and kinase activities.

## REFERENCES

1. Takai, Y., Kishimoto, A., Iwasa, Y., Kawahara, Y., Mori, T. and Nishizuka, Y. 1979. Calcium-dependent activation of a multifunctional protein kinase by membrane phospholipids. *J. Biol. Chem.* 254: 3692-3695.
2. Castagna, M., Takai, Y., Kalbuchi, K., Sano, K., Kikkawa, U. and Nishizuka, Y. 1982. Direct activation of calcium-activated, phospholipid-dependent protein kinase by tumor-promoting phorbol esters. *J. Biol. Chem.* 258: 7847-7851.
3. Kikkawa, U., Takai, Y., Tanaka, Y., Miyake, R. and Nishizuka, Y. 1983. Protein kinase C as a possible receptor of tumor-promoting phorbol esters. *J. Biol. Chem.* 258: 11442-11445.
4. Nishizuka, Y. 1984. The role of protein kinase C in cell surface signal transduction and tumour promotion. *Nature* 308: 693-698.
5. Nishizuka, Y. 1984. Turnover of inositol phospholipids and signal transduction. *Science* 225: 1365-1370.
6. Ohno, S., Akita, Y., Hata, A., Osada, S., Kubo, K., Konno, Y., Akimoto, K., Mizuno, K., Saido, T., Kuroki, T. and Suzuki, K. 1991. Structural and functional diversities of a family of signal transducing protein kinases, protein kinase C family; two distinct classes of PKC, conventional cPKC and novel nPKC. *Adv. Enzyme Regul.* 31: 287-303.
7. Olivier, A.R. and Parker, P.J. 1991. Expression and characterization of protein kinase C- $\delta$ . *Eur. J. Biochem.* 200: 805-810.
8. Osada, S., Mizunou, K., Saido, T.C., Suzuki, K., Kuroki, T. and Ohno, S. 1992. A new member of the protein kinase C family, nPKC $\theta$ , predominantly expressed in skeletal muscle. *Mol. Cell. Biol.* 12: 3930-3938.

## CHROMOSOMAL LOCATION

Genetic locus: PRKCG (human) mapping to 19q13.42.

## PRODUCT

PKC  $\gamma$  (h): 293T Lysate represents a lysate of human PKC  $\gamma$  transfected 293T cells and is provided as 100  $\mu$ g protein in 200  $\mu$ l SDS-PAGE buffer.

## RESEARCH USE

For research use only, not for use in diagnostic procedures.

## APPLICATIONS

PKC  $\gamma$  (h): 293T Lysate is suitable as a Western Blotting positive control for human reactive PKC  $\gamma$  antibodies. Recommended use: 10-20  $\mu$ l per lane.


Control 293T Lysate: sc-117752 is available as a Western Blotting negative control lysate derived from non-transfected 293T cells.

PKC  $\gamma$  (C-4): sc-166385 is recommended as a positive control antibody for Western Blot analysis of enhanced human PKC  $\gamma$  expression in PKC  $\gamma$  transfected 293T cells (starting dilution 1:100, dilution range 1:100-1:1,000).

## RECOMMENDED SUPPORT REAGENTS

To ensure optimal results, the following support reagents are recommended:  
 1) Western Blotting: use m-IgG $\kappa$  BP-HRP: sc-516102 or m-IgG $\kappa$  BP-HRP (Cruz Marker): sc-516102-CM (dilution range: 1:1000-1:10000), Cruz Marker™ Molecular Weight Standards: sc-2035, UltraCruz® Blocking Reagent: sc-516214 and Western Blotting Luminol Reagent: sc-2048.

## DATA

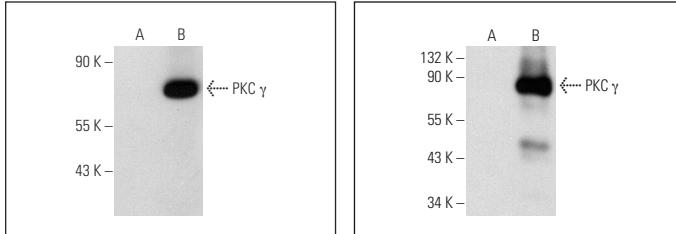

PKC  $\gamma$  (C-4); sc-166385. Western blot analysis of PKC  $\gamma$  expression in non-transfected: sc-117752 (**A**) and human PKC  $\gamma$  transfected: sc-116200 (**B**) 293T whole cell lysates.

PKC  $\gamma$  (D-4); sc-166451. Western blot analysis of PKC  $\gamma$  expression in non-transfected: sc-117752 (**A**) and human PKC  $\gamma$  transfected: sc-116200 (**B**) 293T whole cell lysates.

## STORAGE

Store at -20° C. Repeated freezing and thawing should be minimized. Sample vial should be boiled once prior to use. Non-hazardous. No MSDS required.

## PROTOCOLS

See our web site at [www.scbt.com](http://www.scbt.com) for detailed protocols and support products.