



# SZABO SCANDIC

Part of Europa Biosite

## Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!  
See the following pages for more information!



### Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

### Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

### SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien


T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

[mail@szabo-scandic.com](mailto:mail@szabo-scandic.com)

[www.szabo-scandic.com](http://www.szabo-scandic.com)

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 

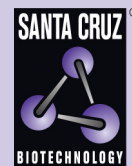

## CUL-2 (m): 293T Lysate: sc-119516

### BACKGROUND

Cullin proteins comprise a distinct family of mediators that participate in the selective targeting of proteins for ubiquitin (Ub)-mediated proteolysis. CUL-1, which is the mammalian homolog of yeast Cdc53, is an integral component of the E3 ubiquitin ligase complex designated SCF. The SCF (Skp1/CUL-1/F-box protein complex) consists of Skp1 associating with both CUL-1 and an F-box protein, such as Skp2, which determines the substrate specificity of the complex. CUL-1-mediated ubiquitination results in the degradation of cell cycle proteins cyclin D, p21 and cyclin E. Another cullin, CUL-3, facilitates the degradation of cyclin E independent of SCF activity, while CUL-2 associates with the tumor suppressing protein VHL and elongin B to form VBC complexes, which structurally resemble the SCF ligase. Proteolysis also occurs by way of CUL-4 associating with Nedd-8, a ubiquitin-like protein, where it too functions as an active component of a multifunctional E3 complex. CUL-5, also designated vasopressin-activated, calcium-mobilizing protein (VACM-1), is also included in the cullin family as it shares substantial sequence homology with CUL-1.

### REFERENCES

1. Kipreos, E.T., et al. 1996. CUL-1 is required for cell cycle exit in *C. elegans* and identifies a novel gene family. *Cell* 85: 829-839.
2. Byrd, P.J., et al. 1997. Identification and analysis of expression of human VACM-1, a cullin gene family member located on chromosome 11q22-23. *Genome Res.* 7: 71-75.
3. Pause, A., et al. 1997. The von Hippel-Lindau tumor-suppressor gene product forms a stable complex with human CUL-2, a member of the Cdc53 family of proteins. *Proc. Natl. Acad. Sci. USA* 94: 2156-2161.
4. Yu, Z.K., et al. 1998. Human CUL-1 associates with the SKP1/SKP2 complex and regulates p21<sup>CIP1/WAF1</sup> and cyclin D proteins. *Proc. Natl. Acad. Sci. USA* 95: 11324-11329.
5. Chen, L.C., et al. 1998. The human homologue for the *Caenorhabditis elegans* CUL-4 gene is amplified and overexpressed in primary breast cancers. *Cancer Res.* 58: 3677-3683.
6. Tyers, M., et al. 1999. One ring to rule a superfamily of E3 ubiquitin ligases. *Science* 284: 601, 603-604.
7. Singer, J.D., et al. 1999. Cullin-3 targets cyclin E for ubiquitination and controls S phase in mammalian cells. *Genes Dev.* 13: 2375-2387.
8. Iwai, K., et al. 1999. Identification of the von Hippel-Lindau tumor-suppressor protein as part of an active E3 ubiquitin ligase complex. *Proc. Natl. Acad. Sci. USA* 96: 12436-12441.
9. Kamura, T., et al. 1999. The Rbx1 subunit of SCF and VHL E3 ubiquitin ligase activates Rub1 modification of cullins Cdc53 and CUL2. *Genes Dev.* 13: 2928-2933.

### CHROMOSOMAL LOCATION

Genetic locus: Cul2 (mouse) mapping to 18 A1.

### PRODUCT

CUL-2 (m): 293T Lysate represents a lysate of mouse CUL-2 transfected 293T cells and is provided as 100 µg protein in 200 µl SDS-PAGE buffer.

### APPLICATIONS

CUL-2 (m): 293T Lysate is suitable as a Western Blotting positive control for mouse reactive CUL-2 antibodies. Recommended use: 10-20 µl per lane.

Control 293T Lysate: sc-117752 is available as a Western Blotting negative control lysate derived from non-transfected 293T cells.

### STORAGE

Store at -20° C. Repeated freezing and thawing should be minimized. Sample vial should be boiled once prior to use. Non-hazardous. No MSDS required.

### RESEARCH USE

For research use only, not for use in diagnostic procedures.

### PROTOCOLS

See our web site at [www.scbt.com](http://www.scbt.com) for detailed protocols and support products.