



# SZABO SCANDIC

Part of Europa Biosite

## Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!  
See the following pages for more information!



### Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

### Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

### SZABO-SCANDIC Handels GmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0


F. +43(0)1 489 3961-7

[mail@szabo-scandic.com](mailto:mail@szabo-scandic.com)

[www.szabo-scandic.com](http://www.szabo-scandic.com)

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic)



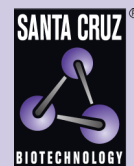

# MAGIX siRNA (m): sc-149228

## BACKGROUND

MAGIX (MAGI family member, X-linked), also known as JM10 or PDZX, is a 334 amino acid protein that contains one PDZ (DHR) domain and exists as three alternatively spliced isoforms. The gene encoding MAGIX maps to human chromosome Xp11.23. Chromosome X consists of about 153 million base pairs and nearly 1,000 genes. The combination of a X and Y chromosome lead to normal male development while two copies of X lead to normal female development. There are a number of conditions related to an unusual number and combination of sex chromosomes being inherited, including Turner's syndrome, Klinefelter's syndrome and Triple X syndrome. Color blindness, hemophilia, and Duchenne muscular dystrophy are well known X chromosome-linked conditions which affect males more frequently as males carry a single X chromosome.

## REFERENCES

1. Givens, J.R., Wilroy, R.S., Summitt, R.L., Andersen, R.N., Wiser, W.L. and Fish, S.A. 1975. Features of Turner's syndrome in women with polycystic ovaries. *Obstet. Gynecol.* 45: 619-624.
2. Bernardino-Sgherri, J., Flagiello, D. and Dutrillaux, B. 2002. Overall DNA methylation and chromatin structure of normal and abnormal X chromosomes. *Cytogenet. Genome Res.* 99: 85-91.
3. Ozcelik, T. 2002. Uncovering the complex mysteries of mosaicism. *Nature* 417: 588.
4. Muntoni, F., Torelli, S. and Ferlini, A. 2003. Dystrophin and mutations: one gene, several proteins, multiple phenotypes. *Lancet Neurol.* 2: 731-740.
5. Deeb, S.S. 2005. The molecular basis of variation in human color vision. *Clin. Genet.* 67: 369-377.
6. Bojesen, A., Kristensen, K., Birkebaek, N.H., Fedder, J., Mosekilde, L., Bennett, P., Laurberg, P., Frystyk, J., Flyvbjerg, A., Christiansen, J.S. and Gravholt, C.H. 2006. The metabolic syndrome is frequent in Klinefelter's syndrome and is associated with abdominal obesity and hypogonadism. *Diabetes Care* 29: 1591-1598.
7. Hayashi, T., Kubo, A., Takeuchi, T., Gekka, T., Goto-Omoto, S. and Kitahara, K. 2006. Novel form of a single X-linked visual pigment gene in a unique dichromatic color-vision defect. *Vis. Neurosci.* 23: 411-417.
8. Augui, S., Fillion, G.J., Huart, S., Nora, E., Guggiari, M., Maresca, M., Stewart, A.F. and Heard, E. 2007. Sensing X chromosome pairs before X inactivation via a novel X-pairing region of the Xic. *Science* 318: 1632-1636.
9. Rolle, U., Linse, B., Glasow, S., Sandig, K.R., Richter, T. and Till, H. 2007. Duodenal atresia in an infant with triple-X syndrome: a new associated malformation in 47,XXX. *Birth Defects Res. Part A Clin. Mol. Teratol.* 79: 612-613.

## CHROMOSOMAL LOCATION

Genetic locus: Magix (mouse) mapping to X A1.1.

## PROTOCOLS

See our web site at [www.scbt.com](http://www.scbt.com) for detailed protocols and support products.

## PRODUCT

MAGIX siRNA (m) is a pool of 3 target-specific 19-25 nt siRNAs designed to knock down gene expression. Each vial contains 3.3 nmol of lyophilized siRNA, sufficient for a 10  $\mu$ M solution once resuspended using protocol below. Suitable for 50-100 transfections. Also see MAGIX shRNA Plasmid (m): sc-149228-SH and MAGIX shRNA (m) Lentiviral Particles: sc-149228-V as alternate gene silencing products.

For independent verification of MAGIX (m) gene silencing results, we also provide the individual siRNA duplex components. Each is available as 3.3 nmol of lyophilized siRNA. These include: sc-149228A, sc-149228B and sc-149228C.

## STORAGE AND RESUSPENSION

Store lyophilized siRNA duplex at -20° C with desiccant. Stable for at least one year from the date of shipment. Once resuspended, store at -20° C, avoid contact with RNAses and repeated freeze thaw cycles.

Resuspend lyophilized siRNA duplex in 330  $\mu$ l of the RNase-free water provided. Resuspension of the siRNA duplex in 330  $\mu$ l of RNase-free water makes a 10  $\mu$ M solution in a 10  $\mu$ M Tris-HCl, pH 8.0, 20 mM NaCl, 1 mM EDTA buffered solution.

## APPLICATIONS

MAGIX siRNA (m) is recommended for the inhibition of MAGIX expression in mouse cells.

## SUPPORT REAGENTS

For optimal siRNA transfection efficiency, Santa Cruz Biotechnology's siRNA Transfection Reagent: sc-29528 (0.3 ml), siRNA Transfection Medium: sc-36868 (20 ml) and siRNA Dilution Buffer: sc-29527 (1.5 ml) are recommended. Control siRNAs or Fluorescein Conjugated Control siRNAs are available as 10  $\mu$ M in 66  $\mu$ l. Each contain a scrambled sequence that will not lead to the specific degradation of any known cellular mRNA. Fluorescein Conjugated Control siRNAs include: sc-36869, sc-44239, sc-44240 and sc-44241. Control siRNAs include: sc-37007, sc-44230, sc-44231, sc-44232, sc-44233, sc-44234, sc-44235, sc-44236, sc-44237 and sc-44238.

## RT-PCR REAGENTS

Semi-quantitative RT-PCR may be performed to monitor MAGIX gene expression knockdown using RT-PCR Primer: MAGIX (m)-PR: sc-149228-PR (20  $\mu$ l). Annealing temperature for the primers should be 55-60° C and the extension temperature should be 68-72° C.

## RESEARCH USE

For research use only, not for use in diagnostic procedures.