



SZABO SCANDIC

Part of Europa Biosite

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!
See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

mail@szabo-scandic.com

www.szabo-scandic.com

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 

SH5-07

Chemical Properties


CAS No. : 1456632-41-9

Formula: C₂₉H₂₈F₅N₃O₅S

Molecular Weight: 625.61

Appearance: no data available

Storage: Powder: -20°C for 3 years | In solvent: -80°C for 1 year

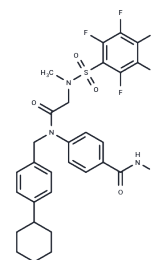

Biological Description

Description	SH5-07 is a hydroxamic acid-based Stat3 inhibitor (IC ₅₀ : 3.9 μM).
Targets(IC ₅₀)	STAT
In vitro	SH5-07, a hydroxamic acid analog of BP-1-102, selectively inhibits Stat3 activity in a dose-dependent manner (IC ₅₀ : 3.9±0.6 μM in vitro), specifically targeting Stat3:Stat3 DNA-binding over Stat1:Stat3, with minimal impact on Stat1:Stat1. By binding to Stat3, it disrupts its interaction with growth factor receptors, thus inhibiting Stat3 phosphorylation. This inhibition leads to reduced expression of Stat3-regulated genes, including Bcl-xL, Bcl-2, c-Myc, Survivin, Cyclin D1, and Mcl-1, following a 24-hour exposure to 5 μM SH5-07, ultimately exerting antitumor effects against cells with constitutively active Stat3.
In vivo	Tail vein injection or oral gavage delivery of SH5-07 inhibits the growth of 90-150 mm ³ established subcutaneous mouse xenografts of human glioma (U251MG) and breast (MDA-MB-231) tumors with aberrantly-active Stat3, associated with decreased Mcl-1, c-Myc, and Cyclin D1 expression. No significant changes in body weights, blood cell counts, organ gross anatomy, or signs of toxicity are observed.
Cell Research	Cells are treated with 0-8 μM agent for 24-48 h. For cell cycle profile analysis, cells are harvested and fixed with 70% ice-cold ethanol and stained with propidium iodide (PI). For apoptosis analysis, cells are collected and stained with FITC-Annexin V using the Apoptosis Detection Kit. Both the DNA content of cells and the Annexin V-positive cells are analyzed by the flow cytometer. Cell cycle phase distribution is analyzed using the Cell-Fit program. Data acquisition is gated to exclude cell doublets.
Animal Research	Mice are injected subcutaneously in the left flank area with U251MG cells in 200 μL of PBS/Matrigel matrix, or MDA-MB-231 cells in 100 μL of PBS. Mice with tumors of 90-150 mm ³ (MDA-MB-231) or 150 mm ³ (U251MG) are grouped for identical mean tumor sizes, administered 3, 5 or 6 mg/kg SH5-07 via oral gavage daily or tail vein injection every 2 or 3 days, and monitored every 3-7 days. Tumor sizes are measured with calipers and converted to tumor volume.

Solubility Information

A DRUG SCREENING EXPERT

Solubility	DMSO: 50 mg/mL (79.92 mM),Sonication is recommended. H2O: Insoluble, (< 1 mg/ml refers to the product slightly soluble or insoluble)
------------	--

Preparing Stock Solutions

	1mg	5mg	10mg
1 mM	1.5984 mL	7.9922 mL	15.9844 mL
5 mM	0.3197 mL	1.5984 mL	3.1969 mL
10 mM	0.1598 mL	0.7992 mL	1.5984 mL
50 mM	0.032 mL	0.1598 mL	0.3197 mL

Please select the appropriate solvent to prepare the stock solution, according to the solubility of the product in different solvents. Please use it as soon as possible.

Reference

Yue P,et al. Hydroxamic Acid and Benzoic Acid-Based STAT3 Inhibitors Suppress Human Glioma and Breast Cancer Phenotypes In Vitro and In Vivo. Cancer Res. 2016 Feb 1;76(3):652-63.

Inhibitor · Natural Compounds · Compound Libraries · Recombinant Proteins

This product is for Research Use Only· Not for Human or Veterinary or Therapeutic Use

Tel:781-999-4286 E_mail:info@targetmol.com Address:36 Washington Street,Wellesley Hills,MA 02481