



# SZABO SCANDIC

Part of Europa Biosite

## Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!  
See the following pages for more information!



### Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

### Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

### SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

[mail@szabo-scandic.com](mailto:mail@szabo-scandic.com)

[www.szabo-scandic.com](http://www.szabo-scandic.com)



[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 



# VectaPlex™ Antibody Removal Kit

## VRK-1000

### Product Images

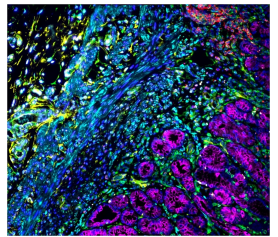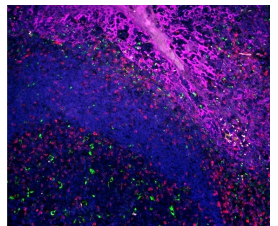

## Short Description

---

VectaPlex Antibody Removal kit is a ready-to-use kit designed to remove antibodies and other non-covalently bound detection reagents from formalin-fixed paraffin-embedded (FFPE) tissue sections in immunofluorescence or tyramide signal amplification staining workflows.

## Additional Information

---

Unit Size	1 kit
Applications	Immunofluorescence, Tyramide Signal Amplification
Recommended Usage	(15–25° C)
Target Species	Universal
Format	Ready-to-Use



# VectaPlex™ Antibody Removal Kit

<b>Cat. No.</b>	VRK-1000
<b>Storage</b>	Ambient (15–25° C)
<b>Description</b>	This ready-to-use kit is designed to remove antibodies and other non-covalently bound detection reagents from formalin-fixed paraffin-embedded (FFPE) tissue sections in immunofluorescence or tyramide signal amplification staining workflows.

## Components

Product Name	Volume
VectaPlex™ Antibody Removal Reagent A	25 ml
VectaPlex™ Antibody Removal Reagent B	25 ml

## Instructions for Use


### Option 1: Immunofluorescence Workflow

#### Preparation for Procedure

Complete standard immunofluorescence staining procedure.  
Mount with a non-hardening (non-setting) mounting medium and image slide.

#### Procedure

1. Remove coverslip from tissue section and wash thoroughly with PBS pH 7.4 to remove residual mounting medium.
2. Add sufficient VectaPlex™ Reagent A to cover tissue section completely (typically ~ 150–200 µL) and incubate for 15 minutes at room temperature.
3. Wash slides briefly with PBS.
4. Add sufficient VectaPlex™ Reagent B to cover tissue section completely (typically ~ 150–200 µL) and incubate for 15 minutes at room temperature.
5. Wash slides with PBS for 5 minutes.
6. The tissue section is now ready to begin another immunofluorescence staining procedure with a protein blocking step followed by the next round of antibody staining, mounting, and imaging.
7. Repeat until the desired number of antigens have been detected.

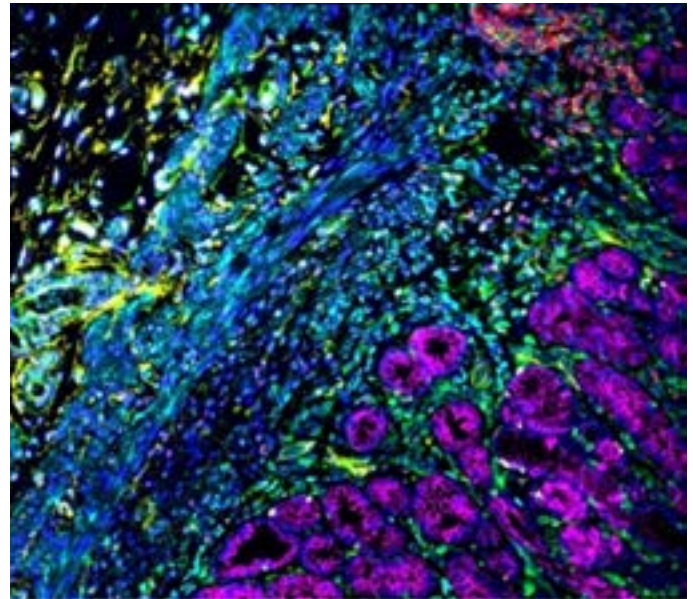

Stomach section stained with 6 mouse primary antibodies (CD20 (Red), CD34 (Yellow), DES (Cyan), AE1/AE3 (Magenta), CD3 (Grey), VIM (Green)) detected with DyLight™ 488 horse anti-mouse IgG and mounted with VECTASHIELD® PLUS with DAPI. This composite image was generated from 6 rounds of immunofluorescent staining using the VectaPlex™ Antibody Removal Kit.

### Option 2: Tyramide Signal Amplification Workflow

#### Preparation for Procedure

Complete standard tyramide signal amplification staining procedure.

#### Procedure

1. Add sufficient VectaPlex™ Reagent A to cover tissue section completely (typically ~ 150–200 µL) and incubate for 15 minutes at room temperature.
2. Wash slides briefly with PBS.
3. Add sufficient VectaPlex™ Reagent B to cover tissue section completely (typically ~ 150–200 µL) and incubate for 15 minutes at room temperature.
4. Wash slides with PBS for 5 minutes.
5. Tissue sections are now ready to begin another fluorescent tyramide signal amplification staining procedure with a protein blocking step followed by the next round of antibody staining using a differently colored dye.
6. Repeat until the desired number of antigens have been detected.
7. Once staining is completed, mount tissue sections and image.

# VectaPlex™ Antibody Removal Kit

## Notes:

1. **Mixing VectaPlex reagents with bleach may result in the formation of a hazardous gas.** For this reason, do not mix reagents or waste solutions from this process with bleach, and use reagents in a well-ventilated area using appropriate laboratory personal protective equipment (PPE).
2. When performing the immunofluorescence workflow, images from successive staining cycles can be registered and overlaid using free software such as Fiji (<https://fiji.sc/>).
3. When performing the immunofluorescence workflow, use of a mounting medium containing DAPI during each staining cycle facilitates image alignment.
4. This product has been developed for use across a wide range of FFPE sample and staining conditions. Up to six rounds of successive staining and antibody removal have been demonstrated. However, performance may vary depending on the sample or experimental design used. For this reason, control experiments using single-plex staining with and without VectaPlex usage are recommended.