



SZABO SCANDIC

Part of Europa Biosite

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!
See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

mail@szabo-scandic.com

www.szabo-scandic.com

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 

SARS-CoV-2 (COVID-19) Envelope FFPE 293T cell pellet Block

Cat No. GTX435642

Application	IHC-P (cell pellet)
Reactivity	SARS Coronavirus 2
Species	SARS Coronavirus 2

Package
1 block

PRODUCT

Summary

Block size is 2.2 x 2.2 x 0.5 cm (expected to yield ~40-50 sections of ~4-5 μ m thickness), but can vary slightly. Cell pellet consists of 293T cells overexpressing SARS-CoV-2 Envelope (C-terminal 3x DDDDK-tagged) protein fixed with 4% PFA (15 min, room temp) 48 hrs post-transfection, and then embedded in paraffin.

APPLICATION

Application Note

*Optimal dilutions/concentrations should be determined by the researcher.

Suggested dilution	Dilution
IHC-P (cell pellet)	Assay dependent

Not tested in other applications.

PROPERTIES


Storage

Store at room temperature.

Note

For laboratory use only. Not for any clinical, therapeutic, or diagnostic use in humans or animals. Not for animal or human consumption.

DATA IMAGES

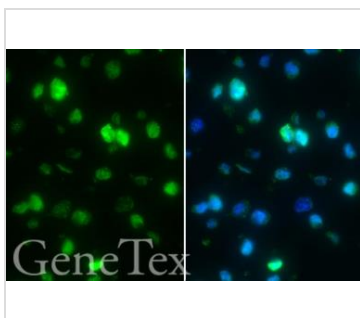

GTX435642 IHC-P (cell pellet) Image

DDDDK tag antibody detects DDDDK-tagged SARS-CoV-2 (COVID-19) envelope protein by immunohistochemical analysis.

Sample: Paraffin-embedded SARS-CoV-2 (COVID-19) Envelope FFPE Cell Pellet Block (GTX435642).

Green: DDDDK-tagged SARS-CoV-2 (COVID-19) envelope protein stained by DDDDK tag antibody (GTX115043) diluted at 1:500.

Blue: Fluoroshield with DAPI (GTX30920).

Antigen Retrieval: Citrate buffer, pH 6.0, 15 min



For full product information, images and publications, please visit our [website](https://www.genetex.com).